

2023 CATALOGUE



Welcome to the **Equi-Still Stocks and Dental Instrumentation Catalogue!**

At Equi-Still, we sell safety. Whether it's stocks or high quality equine dental instruments, we've got your back!

Stocks:

For dentistry, reproduction, radiology, or standing surgery, our *portable or stationary stocks*' three-pole, hot-dipped galvanized steel and aluminium construction, pressure-independent quick-release latching mechanism, plus removable/washable padding in crucial areas keep horses, vets, and handlers safe. Designed for safety, robustness and ease of use with many options.

Instruments:

To complement NZVA's Equine Dentistry courses, we are importing a wide range of instrumentation, with *exquisite quality instruments* for precision examination, floating, extraction, diastema therapy and more, for your world class dentistry! From *HDE*, *Pegasos*, *MAI*, and others!

We also offer a line of budget essential instruments which will allow you to provide basic diagnostic and therapeutic dental care while proving the viability of dental practice expansion... before investing in more pricey advanced dental instrumentation solutions. For more details, see page 27.

Let us help you today!



Contents

Page 4 Product Comparisons
Page 9 Equi–Still Portable Stocks
Page 16 Equi–Still Instruments
Page 33 Pegasos of Germany
Page 53 HDE of France
Page 106 MAI of USA

Product Comparisons

FLUSHING INSTRUMENTS

Rinsing Syringe

400 ml

Flushing pistol, for one-handed use 200 ml

A CONTRACTOR OF THE PARTY OF TH

HDE

The HydroJet 2 (DISCONTINUED)



Pegasos Equus Line

The Equus Jet Ecowith lead battery

- High pressure (up to 11 bars)
- Flow regulated via the wands
- Lead battery (up to a week autonomy)
- Fastconnecting and and and Rinsing Wand
- Double filter system (external and internal)

The Equus Jet the premium unit

As for the Equus Jet Eco, but with additional features:

- 5 set levels of pressure regulated via the white control panel on top
- Lithium ion battery (up to one month autonomy)
- Optional extra: wired foot pedal for flow / pressure regulation



MAI

The Edge Equine Water Pick Kits

- Quick ChangeWater Pick & FlushingWand Kit with Bottle or
- Quick Change Water Pick & Flushing Wand with Tubing Assembly to fit on hose
- Great for flushing the mouth, but doesn't achieve the pressure of Equus Jet







DENTAL SCOPES

Intraoral Camera

This one is more fragile, but will let people get into a camera for a very affordable price.

The camera is on page 15 of the catalogue.







DentalscopeCamera by HDE

Exquisite new camera. Please see catalogue page 101

Starter Kit ACC_CD

Starter Kit +Wireless USB
Connector

ACC_CDUSB

Starter Kit +Wireless monitor ACC_CDM





MAI SPECULUMS

X Spec[™]

Next generation equine speculum, an enhancement of the original AlumiSpec $^{\text{TM}}$, known and trusted throughout the industry.

- Aircraft grade aluminum with a durable powder coat finish
- Wider design to accommodate digital radiographic plates
- Recommended for use with sedated animals
- Tested on hundreds of patients in an equine dental practice
- Net weight of 2.5 kg
- Replaceable synthetic straps
- Patented design

AlumiSpec™

- Replaceable straps in "cinch-lock" locking system allows infinite degrees of access to the oral cavity
- Easy to clean and is pliable at all temperatures
- Lightweight aluminum, net weight of 2.8 kg
- Durability not compromised by lightweight design



MAI SPECULUM REPLACEMENT PARTS

53602 ReplacementOpeningLatch Strap for AlumiSpec™ & X Spec™

53603 Cam Buckle for Opening Latch Strap for AlumiSpec™ & X

Speco® ReplacementX Spec™ Bite Plate Pads with Holes (Set of 2)

53676 ReplacementPoll Strap for AlumiSpec™ & X Spec™

53675 Bolt Kit for AlumiSpec[™] & X Spec[™] –for One Side

53425–1Magnetic Speculum Plate (for use if you use a magnetic spec light, as magnetswon't stick to aluminium.

HDE Speculums

2000 series Millenium Speculum with mouth plates and Biothane straps

SpeculumOPTIMUM (without plates)



Speculum Standard (without plates)
Pony Speculum Standard (without plates)



SpreaderarchforOPTIMUM & MILLENIUM Holds speculum arms wide. wonderful.



Cheek Retractor



Pony Speculum plates (pair)

Speculum plates (all are pair) Classic plates standard



Incisor plates, rubberised




Incisor malocclusion plates standard horse



For the Optimum, I prefer theseones:

Incisor malocclusion plates

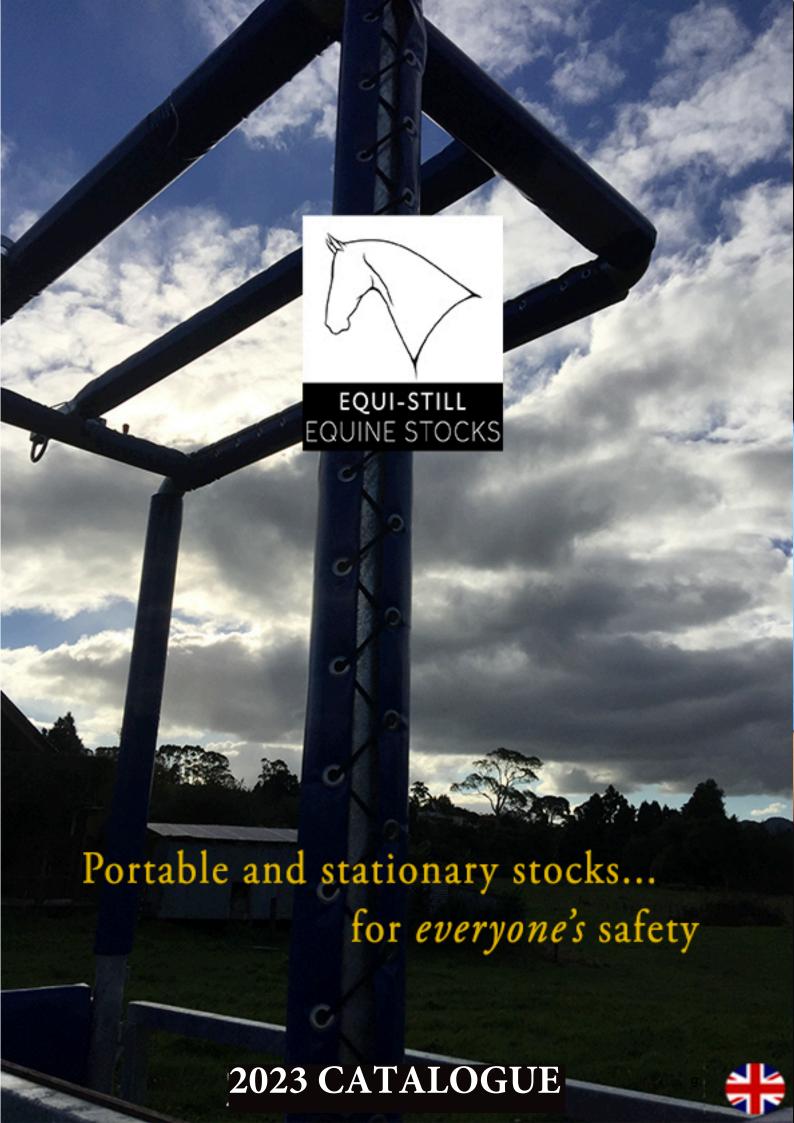



PORTABLE STOCKS
TUI
KEA

STATIONARY STOCKS
TUI
KEA

OPTIONAL EXTRAS

MINI-BAR
EQUI-TABLE
EQUI-TIE
EQUI-COVER

COMMERCIAL GUARANTEE



#1 New Zealand manufacturer of equine portable and stationary stocks

PORTABLE STOCKS

designed for safely, robustness and ease of use by a U.C. Davis-trained equine veterinarian and ICU Technician heavy, hot-dipped galvanized steel

best quality closed cell foam padding with heavy duty woven polyester coated on both sides with flexible PVC babsence of a fourth post

>quick-release rear latch

winch, water tank, and pump, powered via a jump-start battery pack, allowing you to use it with multiple practice vehicles without specialty wiring

➤ one *Equi-Table* is included, which may be mounted on one of three corners

ristocks tip up and down via winch, and may be raised and lowered even when not attached to vehicle

>patients may be restrained in forward or rearward configuration to assist in dental, reproductive, radiological, or surgical procedures

best quality engineering: which means our stocks back and tow beautifully, barely tracking around bends our make and models are registered with the USA NHTSA (National Highway Traffic Safety Authority) stocks met voluntary NATM Compliance Certification in 2018, meaning they met or exceeded most USA

a lasting investment in your practice and your safety

MODEL TUI

• one side of stocks is formed from square steel bars, allowing better access to patient

states' requirements for registration

- The other side is heavy plywood, which may be drilled to allow rope placement to fit your smaller patients
- a rope, or Optional Extra *Equi-Tie*, holds patient in place, whether horse is placed head-first or rear-first for reproduction work









PORTABLE STOCKS

MODEL KEA

- two sides of stocks are solid steel with holes for bars
- a rope, or Optional Extra *Equi-Tie*, may be used instead to hold patient in place, whether horse is placed into stocks head-first, or rear-first
- the holes in the steel sides allow use of bars to hold horse in place, should you desire them









STATIONARY STOCKS

- designed for safely, robustness and ease of use by a U.C. Davis-trained equine veterinarian and ICU Technician
- ➤ heavy, hot dipped galvanized steel
- best quality closed cell foam padding with heavy duty woven polyester coated on both sides with flexible PVC
- absence of a fourth post
- > quick-release rear latch, which combined with only three posts, means stocks open up fully in emergency
- > stationary stocks come with separate mounting pipes to inset into floor concrete where you wish the stocks to be placed, so the stocks may later be moved
- > stocks include one *Equi-Table*
- > stationary stocks have no floor, winch, pump, water tank, etc.
- > a lasting investment in your practice and your safety

TUI STATIONARY

- as with our portable stocks, one side of stocks is formed from square steel bars, allowing better access to patient
- a rope, or Optional Extra *Equi-Tie* holds the patient in place, whether horse is placed head-first or rear-first for reproduction work
- best quality closed cell foam padding with heavy duty woven polyester coated with flexible PVC on both sides
- absence of a fourth post
- quick-release rear latch
- a lasting investment in your practice and your safety



KFA STATIONARY

- two sides of stocks are solid steel with holes for bars
- a rope, or Optional Extra *Equi-Tie*, may be used instead to hold patient in place, whether horse is placed into stocks head-first, or rear-first
- the holes in the steel sides allow use of bars to hold horse in place, should you desire them



OPTIONAL EXTRAS

THE MINI BAR

- mounted on the lowest side rail of the stocks to help contain your smaller patients, be they equine, camelid, bovine or anything else
- · easy on, easy off






THE EQUI-TIE

- heavy, thick, soft nylon rope with brass snap
- keeps patients in place during procedures
- plywood backing of *Tui* model may be drilled at heights appropriate to your patients' sizes to run the soft rope through

THE EQUI-TABLE

- removable table, one provided with stocks
- additional tables may be purchased, to your specs
- mounts on any of three different corners of stocks
- designed to tip to read ultrasound or other screens

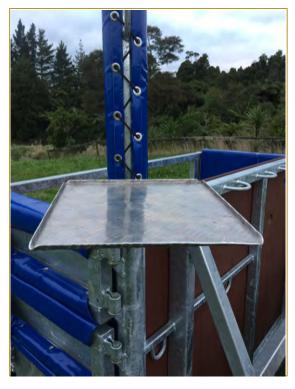



THE EQUI-COVER

- stocks cover for use while on the road or for storage
- same poly/PVC as stocks upholstery

OUR COMMERCIAL GUARANTEE

1. SCOPE

In additiontotheLegalGuarantee,EQUI-STILL PORTABLE EQUINE STOCKS (EQUI-STILL) offers practitioners an EQUI-STILL Commercial Guaranteewhich applies to all new EQUI-STILL-branded portable and stationary stocks purchased by the practitioner from EQUI-STILL, from Aug 20, 2020.

2. DEFINITION

a) Product Commercial Guarantee

EQUI-STILL guarantees to the practitioner against breakagehis/herEQUI-STILL-branded portable and stationary stocks for a 12-month period following its date of purchase from EQUI-STILL. This applies to:

- -Broken parts (except "wear and tear" parts paragraph 3.b)
- -Labour costs
- -Transport costs.

b) After-SalesServiceCommercial Guarantee

All repairs made by EQUI-STILL on EQUI-STILL portable and stationary stocks are guaranteed for 6 months. This guarantee covers:

- Broken parts replaced in the last 6 months (except "wear and tear" parts paragraph 3.b)
- Labour and transport costs for the replacement of these parts

c) Free annual check-up

After the initial 12-month period of guarantee, EQUI-STILL and our Official Technical Centers* offer one Free Annual Check-up to extend the life of your equipment. WOF requirements in New Zealand are the responsibility of the purchaser. The Check-up includes the following:

- -Complete test of the stocks
- -Cleaning
- -Greasing

Labour costs are free of charge, transport costs are at the practitioner's expense.

This free check-up is valid once once every calendar year.

Further to this check-up,in case repairs are required, spare parts, repair labour and transportation costs are at the practitioner's expense. In this case, an estimate will be sent for approval before any repair is done.

3. IMPLEMENTATION

a) Product registration

The product commercial guarantee, the after-sales service commercial guarantee and the annual check-upwill only be implemented if EQUI-STILL products are registered on EQUI-STILL's website within 15 days following the purchaseby the practitioner from EQUI-STILL or one of its official distributors.

b) Disclamer

The commercial guarantee applies to all EQUI-STILL Portable Equine Stocks except:

A - The following parts:

- -tyres, (other than manufacturer's warranty on defects).
- jump-startbattery (other than manufacturerr's warranty) .
- electrical system or electrical components, if any alteration to the electrical power system or power items have been performed other than by an Official EQUI-STILL Service Center. This includes, but is not limited to, power input source, water pump, winch, lights, or any wiring.
- B If the problem is caused by an accident
- C If the damage is caused by the user, including driving at speed over irregular roads or speed bumps.
- D If the user fails to comply with the instruction manual and/or EQUI-STILL cleaning and service recommendations
- E If the VIN number of the product is removed or modified
- F If the productis modified by anyone other than an Official EQUI-STILL Service Center or if it is modified with non-EQUI-STILL parts.
- G If the productwas not registered by the practitionar on EQUI-STILL's website within 15 days after its purchase



Making dentistry fun again!







DENTAL HALTER

DENTAL HEAD SUPPORT

INTRAORAL CAMERA

HEAD LIGHT SET

HAND INSTRUMENTS

LONG WOLF TOOTHELEVATORS
DENTAL MIRRORS

K-FILES



EQUI-STILL: New Zealand's #1 manufacturer of equine portable and stationary stocks



>Round Mirror on T-Handle

•smaller mirror may be easier to handle in the mouth

➤ Mirror Spares

•to replace mirror in case of breakage

>Alligator Forceps

•40 mm long, with 2-3 cm tips at 90°.

➤ Double Ended Scaler

•perfect for scaling tartar from incisors and canines, plus evacuating rostral diastemata

➤ Periodontal Probe with Black Increment Markings

•for determining and documenting depth of periodontal pockets

>Large Curette Scaler

•sturdy instrument for scaling teeth, evacuating diastemata, and as a general pick

➤ Dental Pick

•sturdy instrument for evacuating diastemata and as a general pick—no separate pieces to lose

>Explorer Fine Tipped

•Fine tipped probe for investigating for pulp exposures and other pathologies—remember, this one is more fragile than the pick and the scaler, above

➤ Leather-Look Synthetic Zippered Case

•professional looking, with easy clean exterior



Finally, a safe and convenient support for dental and other treatments of the equine head!

Made from heavy duty Biothane, which remains soft and pliable, no matter how cold the temperature, no matter how many washes and disinfections you give it!

➤ Heavy duty snaps, rings, and buckles are made from high quality stainless steel

All parts are replaceable, for future-proofing your dental halter

➤ Chicago Screws (rather than rivets) throughout allow easy removal or replacement of straps

Designed for stable support of head--no tipping from side to side if horse moves head off centre--on the head rest, so no need to continually re-position halter and patient

➤ Snap hook on chin strap for quick on and off fitting, no need to unbuckle the crownpiece

Soft cushion on head rest is coated in rubber for easy cleaning. This head rest may easily be replaced as it wears over time, whereas the halter itself has a long life span. You will get a discount on the new replacement head rest when you return your old one!







Our own standard head support lets you raise or lower your patient's head with a simple one-handed cleat.

- ➤ The Equi-Support system includes 4 metres of soft, heavy nylon braided rope and a quick-release cleat for safety
- Easy, one-handed action lets you release the support as quickly or slowly as you wish
- > Dental halter not included





DENTAL HEAD SUPPORT



Our dental head support suspension system allows for a simple, smooth way of raising or lowering the horses head to the exact height needed.

- The support system includes strong stainless steel snaps and quick release latch for safely
- The quick release rope should be positioned to the side of the working space (e.g. stocks, doorway), allowing the practitioner to step away from the horse while being able to instantly release the head support when needed





INTRAORAL CAMERA-SCOPE

This brilliant little number lets you photograph and video on your smartphone!

- ➤ Its articulating head bends more than 90 degrees, enabling you to inspect deep in the alveolus
- > Stainless steel shaft, with an easy-clean, comfortable handle.
- Self-containedWiFi systemlets you connect multiple nearby analogue or iOS phonesor iPads for viewing by you and clients—capture images or videos using your device's photo capability
- > Excellent detail and colour of image
- > Six variable brightnessLEDs let you see!
















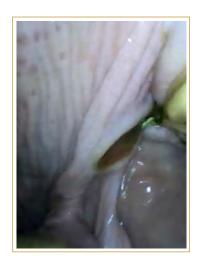





NSE USB VIDEOSCOPE



MODEL UE-100 HAS BEEN DESIGNED SPECIFICALLY FOR EQUINE DENTISTRY IN CONSULTATION WITH A TOP EQUINE DENTAL VET.!

➤ Just plug in and play for real-time observation and recording of your oral (or other) endoscopic exam

Multi functions: Photo, Video, Freeze. Zoom in, zoom out.,

➤ Review image 180° reversible.

➤ Wand with camera is 360° rotatable

1/9" COMS Sensor,1280*720Video Output1280x720, 30PWorking Length42cmAngle of View70 degreeUSB Cable2.8meter length, USB 2.0

Diameter 8mm (working pole)

Operation Software MedCam










NSE EXPERT HEADLAMP



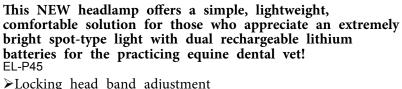

► Locking head band adjustment

➤ Two Li-ion polymer batteries of 2000 mAH power capacity battery 3.7 V, weight ONLY 150 g. (Incl battery)

Easy on-off five-level dimming buttons and battery power indicator lights. Spare battery.

Maximum illumination 40,000 Lux at 200mm. Light spot diameter 88 mm at 450mm





NSE SURGICAL HEADLAMP



This NEW headlamp offers a simple, lightweight, comfortable solution for those who appreciate an extremely bright spot-type light with dual rechargeable lithium batteries for the practicing equine dental vet! EL-S50100

► Locking head band adjustment

➤ Two Li-ion polymer batteries of 2000 mAH power capacity battery 3.7 V, weight ONLY 150 g. (Incl battery)

Easy on-off five-level dimming buttons and battery power indicator lights. Spare battery.

Maximum illumination 50,000 Lux at 200mm. Light spot diameter 100 mm at 500mm





HEADLIGHT SET



This headlamp has many features designed to address issues often found with headlamp in a high-comfort setup for the practicing equine vet.

Locking ratchet head band adjustment system.

➤ PATENT PENDING battery holder evenly distribute weight of the battery, balancing the light and preventing slipping of the head lamp when the practitioner tilts their head or sweats.

➤3400 mAH power capacity battery (average batteries 2600mAH) with power level indicator and rotary dimmer for your own preferred light intensity

➤ Heavy-duty on-off toggle switch


Long light housing armature leaves adequate room for safety glasses and to helps minimise lens fogging







- ➤ 3 pivot-point screws let you adjust light beam exactly where you want it, without repeated adjustments
- Lens reflects light evenly both across the beam width, and from the front and the back of the mouth as well
- ➤ Up to 1000 lumens clear WHITE light with no tinting, offering accurate tissue coloration
- ∠Back up light (for that day you leave your battery charger at home)
- Dipped rubber coating makes your head lamp water and dust resistant

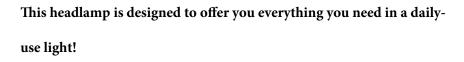






HAND INSTRUMENTS

IONG SHAFTED ELEVATORS

- Solid stainless steel, not hollow
- Nice for wolf tooth extraction

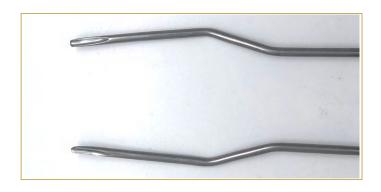



DENTAL MIRRORS

- ➤ Long handled stainless steel mirrors
- > Traditional round and rectangular
- ➤ Each mirror comes with a thick protective sponge cover









MULTIPLE SIZES AVAILABLE

We stock packets of six:

- ➤ single length/ single size
- > multiple sizes/ single length

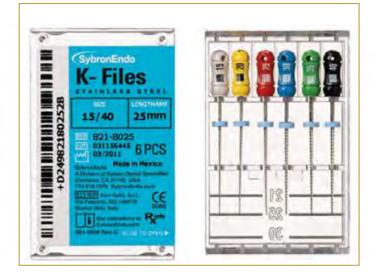




Image courtesy of Mani, Inc.

Size

	15-40	20	30	35	45-80	50
21 mm	A	В				
25 mm			C		D	\mathbf{E}
31 mm				F	G	

MULTIPLE SIZES AVAILABLE

We stock packets of six:

- ➤ single length/ single size
- > multiple sizes/ single length

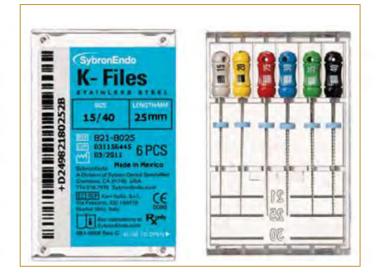




Image courtesy of Mani, Inc.

Size

	15-40	20	30	35	45-80	50
21 mm	A	В				
25 mm			C		D	\mathbf{E}
31 mm				F	G	

Length

BUDGET DIAGNOSTIC AND TREATMENT ESSENTIALS

FOR GETTING STARTED IN EQUINE DENTISTRY

Pony 4 Prong Extraction Forceps Round Handle



Pony 4 Prong Extraction Forceps Ergonomic Handle



Pony Spreader Simple Joint Round Handle



Pony Capps Spec, Biothane Straps with Two Mouth Plates





Pony Millennium Speculum, with Pony standard Mouth Plates & Biothane Straps



Endodontic Plates Half open, set of left and right





Gunther Speculum Excellent for displaced mandible radiographs



Fragment Forceps nose 4 cm **Fragment Forceps** Nose 6 cm Frag Forceps Nose 8cm Slab Fracture Fragment Forceps SET 30 &60 mm Reynolds Cap **Extractor Lower** Reynolds Cap **Extractor Upper** Hewson's Cap Extractor w spring 16"

Incisor Spreader 7'



Incisor forceps 11"



Incisor Fragment Forceps 19"





Instrument Box



Dental Mallet



Slotted Hammer



Alligator Forceps 30cm, Tip length 2cm



Perio Probe Black Bands



Perio Probe Coloured Bands



Dental Mirror Round, Angled



Scaler fine: (Cleaning Probe Short Tip with T- Handle)	
Scaler: Double Point	
Explorer: Probe with Point	
Dental Pick: Probe with Fibre Handle	
Dental Flushing Syringe Plastic 400ml	
Dental Flushing Syringe Brass	

Cheek Retractor Wolf Toothor Root Tip Forceps Slim 31cm Wolf Toothor Root Tip Forceps Slim 20cm



Pegasos Edition

2019 CATALOGUE


High Precision Extraction


the frey the

For veterinary use only



GENERAL EQUIPMENT

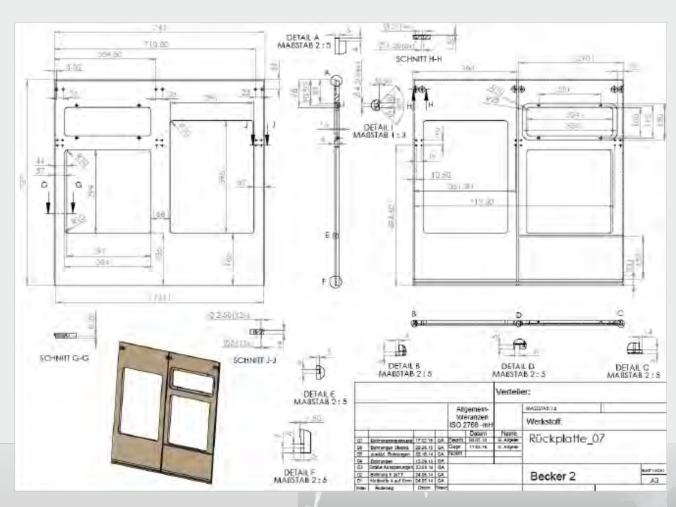

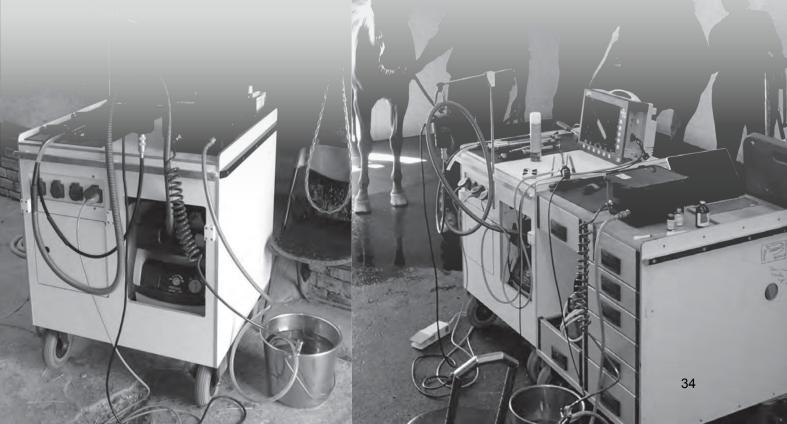

EQUODENT - DENTAL TROLLEY

Our treatment trolley was developed to place the expanding equipment for routine equine dentistry in a safe and tidy way. It creates a workplace which is easy and fast to set up, take down and accommodate in the practice car.

Its versatile cube shape and mixture of materials (coated multiplex boards and

aluminum drawers) make the trolley light and robust at the same time. Its surfaces are durable and easy to clean. It fits through every door with standard widths. It is easy to manage and fits into most compact vans.





ArtNo.#	Description
PEGA_313000	Treatment Trolley
PEGA_210122	Equodent Dental Burr Model Richard Miller - Watercooling, HDE connection
PEGA_221001	Equodent Dental Burr - Watercooling, key drive
PEGA_221002	Equodent Dental Burr -Watercooling, HDE connection
ART- 00015	Waterbottle and connection
PEGA_315050	Diamond coated burr Bit (3,175 mm / 38 mm)
PEGA_220020	Carbide burr (3,175 mm / 38 mm)
PEGA_220021	Carbide burr (3,175 mm / 51 mm)
PEGA_220022	Carbide burr (3,175 mm / 64 mm)
PEGA_220023	Carbide burr (3,175 mm / 77 mm)
PEGA_220000	Storage block for burr bits (3,175 mm), autoclavable (metal, black)

EQUODENT DENTAL BURR MODEL RICHARD MILLER WITH WATERCOOLING

This bullet nose head dental burr is especially useful for horses with smaller heads and for the last cheek teeth (11) where most conventional equipment quickly faces anatomical limitations.



EQUODENT DENTAL BURR WITH WATERCOOLING

This burr with IC 300-Nakanishi head was developed for following applications:

- diastema work
- cheek tooth segmentation
- tooth crown remodeling to improve the grip of the extraction forceps
- pre-drilling during Minimally Invasive Transbuccal Extraction

The shaft bearing is greased which makes it "water-proof" and running quiet, without vibrations.

It has an external water supply enabling watercooling during drilling. This is necessary because the temperatures of the diamond burrcan become very high and can damage the tooth. It is also available for both HDE shaft and dremel Key Drive.

This diamond coated diastema burr with a diameter of 3,175 mm /0.12 inches was developed for diastemas work. It fits the IC-300 Nakanishi head.







CARBIDE BURRS

For segmentation we offer carbide burrs in 4 different length.

For a clear storage, one can use our autoclavable storage block.

DENTAL PROBES



This dental probe was developed for exploring teeth and cleaning interdental spaces. It is available in a long and robust version (18,5 mm / 0.71 inches) and in a short and fine version (15 mm / 0.59 inches).

DENTAL HALTER

The separation of dental halter and mouth gag reduces pressure on the incisors during treatment. Switching between working with and without mouthgag becomes easy. Furthermore it provides a safe handling of challenging horses during injection. The dental halter is covered with leather

The dental halter is covered with leather and has a soft and durable padding.



DENTAL X-RAY BLOCK



This block is a simple tool to open the mouth in 3 different angles for mouth radiographs. It is well tolerated in sligthly sedated horses.

ArtNo.#	Description
PEGA_311001	Dental Probe long strong (18,5 mm)
PEGA_311003	Dental Probe short fine (15 mm)
PEGA_311005	Dental Halter
PEGA_311008	X Ray Block

EQUUS JET RANGE

Multi-purpose rinsing devices

The Equus Jet range has been developed as an autonomous, robust, compact cleaning and cooling system with the mobile equine dentistry practitioners and clinics in mind. It is ideal for cleaning the mouth before and after dental work, diastema, rinsing after tooth extraction, sinus flushing or cooling when working with handpieces.

EQUUS JET

- High pressure (up to 11 bars)
- 5 set levels of pressure regulated via the white control panel on top
- Lithium ion battery (up to one month autonomy)
- Possibility to use a wired foot pedal to gradually regulate flow and pressure
- Fast connecting and disconnecting hoses and wands
- Double filter system (external and internal)







EQUUS JET ECO

- High pressure (up to 11 bars)
- Flow regulated via the wands
- Lead battery (up to a week autonomy)
- Fast connecting and disconnecting hoses and wands
- Double filter system (external and internal)



The Equus jet eco and the Equus jet include the following accessories:

• One flushing wand • One diastema wand • An input polyurethane hose pipe with a strainer and filter • An output spiral polyurethane hose pipe with fast connector on one side and fast connector with valve on the other, allowing to change the wands without shutting off the pump to save time • A special lead or lithium ion battery charger

ArtNo.#	Description
PEGA_EJ-EJE02	EquusJet Eco
PEGA_EJ-CBEP01	ChargerEquus Jet Eco
PEGA_EJ-LDEJ1	Diastema handpiece
PEGA_EJ-LDEJ2	Rinsing handpiece
PEGA_EJH-EJ11	EquusJet
PEGA_EJ-EJ04	Equus Jet - Option Pegasos Trolley
PEGA_EJ-CBLI14	ChargerEquus Jet
PEGA_EJ-LDEJ1	Diastema handpiece
PEGA_EJ-LDEJ2	Rinsing handpiece
PEGA_EJH-PDC02	Equus Jet Foot Pedal
PEGA_EJH-KLS02	Equus Jet Sinus Kit
PEGA_EJ-EJV21 PEGA_FAB-EDH50630 & PEGA_FAB-DH5190	Equus Jet Connect including cooling and peristaltic pumps
PEGA EL-PC01	Equus let Connect Wireless Foot Pedal

INSTRUMENTS FOR DENTAL EXTRACTIONS





THE EXTRACTION SYSTEM FROM PEGASOS4D MODEL: HAUPTNER

Edition: Johann Heinrich Friedrich Günther







Comments on Model and Edition

Model Hauptner: The development and construction of these instruments was supported by the company Hauptner. Hauptner has more than 150 years experience on making instruments and produced in the past some of the best forceps available at that time.

Edition: Johann Heinrich Friedrich Günther


Who was J. H. Friedrich Günther and what influenced id he have to the method of oral extraction of teeth?

He was born in 1794 in Germany. He got educated in Jena, Berlin and Hanover. After he had finished "Royal Equine Vetschool" Hanover in 1818 he started to work as a veterinary surgeon. Only one year later he was called by the school to become a teacher. Apart from his work at the vetschool he ran a private clinic where he tought the older students the practical part. Additionally he was an inspector at the Marstall and gave lessons at the military academy. Since 1840 he did research as well.

In 1847 he became the director of the Hanovarian vetschool. He strongly increased the level of education of veterinary surgeons during that time.

Alongside teaching he did a tremendous high quality work concerning the development and construction of surgical instruments. Furthermore, he gathered a comprehensive collection of diseased teeth. He is, indeed, considered as a founder of the equine dental surgery. Together with his son he wrote a booklet about horses' teeth which was published in 1859, one year after he had died.

Even 50 years after his death his extraction forceps were widely used and highly valued. He was doing oral extraction of equine teeth at a high level. Unfortunately this knowledge got lost over the following century. His early work is only little known and it is a mission to make it accessible again to veterinarians. His early literature of 1859 is highly recommended because it is still up to date in many respects. In order to support the rediscovery of his work we named this forceps and spreader line after him: Edition "Günther".

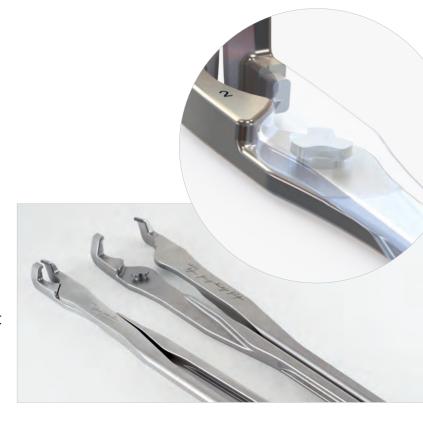


The Spreader and Forceps System was developed to improve oral extraction techniques in equine dentistry.

Great efforts have been invested to design the instruments as ergonomic and efficient as possible. The different sizes of spreader blades enable a standardized stepwise working process. Thereby, the duration and success rate of the extractions are influenced very positively.

The extreme precision of the joint creates the possibility to disassemble and reassemble the single halves of the instruments.

It enables a variety of combinations that has not been available yet, like the combination of spreader and forceps within one instrument. It is most helpful in cases of slab fractures and very tight interdental spaces especially in the lower jaw. It also enables a thorough cleaning of all parts.

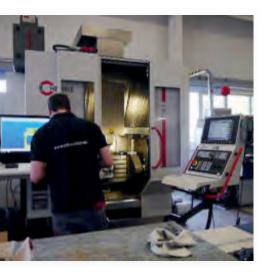

Manufacturing


- Instruments are made of heat-treated surgical stainless steel.
- The body is milled out of one piece on a 5-axis-CNC-milling centre
- The spreader inserts are laser-welded and labeling is done with a laser marking system.
- The construction was tested with FEM-analysis and does still have, at a 80 kg handle pressure, a safety factor of 3.
- Made in Germany*.
- 5-year warranty.

General notes on maintenance

The possibility to disassemble the instruments enables easy and thorough cleaning. This should be done after each use. From time to time the fulcrum pin should be treated with surgical instrument silicone spray.


All instruments are fully autoclavable and chemically disinfectable.

MOLAR SPREADERS

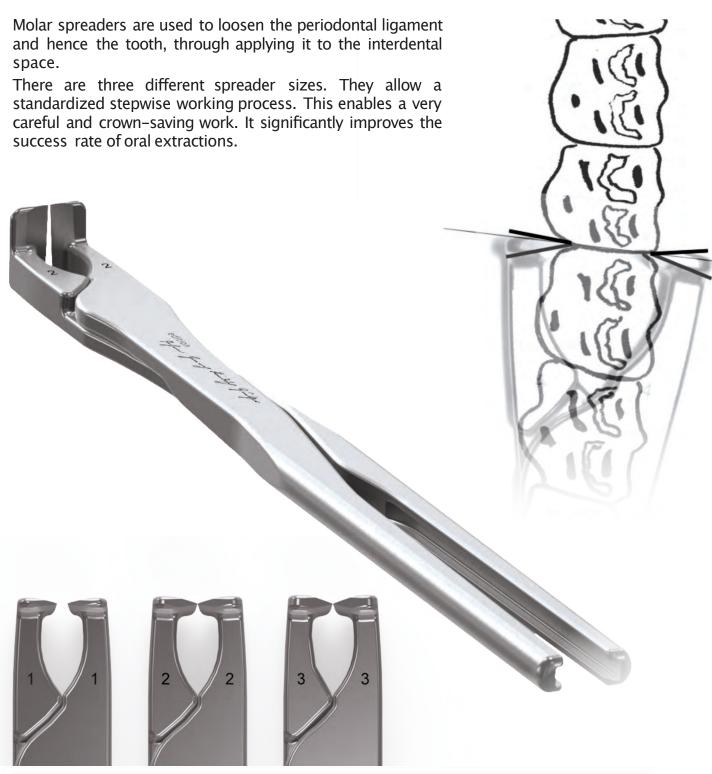

ArtNo.#	Description	
PEGA_320001	Molar Spreaders No. 1 - thickness 4 mm	
PEGA_320002	Molar Spreaders No. 2 - thickness 5 mm	
PEGA_320003	Molar Spreaders No. 3 – thickness 6 mm	

Differences between No. 1–3

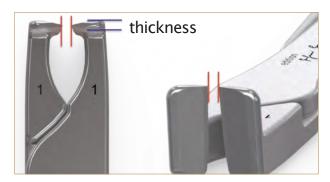

Molar Spreader No.1, the thinnest one with a 4 mm (0.16 inches) thickness (blue), does additionally have a cleft (red) between the jaws. Therefore, the handles don't need to open that far. This enables a more ergonomic first spreading step. One can apply a much more controlled force. This is especially helpful in case of fragile teeth.

The varying jaw thickness enables different spreading effects.

Spreader No. 1 – Narrow – the 4 mm (0.16 inches) thickness works much finer than No. 2 – Medium. No. 1 is mainly used as the first step spreader for the maxillary cheek teeth because they are wider than mandibular cheek teeth.

Spreader No. 2 – Medium – the 5 mm (0.19 inches) thickness works for the upper jaw as a second step and for mandibular cheek teeth as a first step.

Spreader No. 3 – Wide – the 6 mm $(0.24 \, \text{inches})$ thickness is used only partially (1/3-2/3) for the maxillary cheek teeth and up to full effect for mandibular cheek teeth. It is recommended to use the forceps to loosen the tooth with rotational movements between the spreading steps, otherwise the risk of breaking the tooth/rootswill increase.

Interchangeability







Additionally, the application of spreaders in combination with forceps is possible thanks to its interchangeability. This is especially helpful in cases where the interdental space is not wide enough to use conventional spreaders. This is the case in slab fractures and for the lower iaw.

Cavity: one has to be careful on side of the forceps part otherwise breakage of the tooth can occur.

EXTRACTION FORCEPS



Extraction forceps are used to further loosen teeth with rotational and tilting movements and finally to extract the teeth in combination with a fulcrum.

The Extraction Forceps No. 1 is particularly helpful for mandibular cheek teeth extractions but can also be used for the upper jaw.


The Extraction Forceps No. 2 is specifically designed for the wider maxillary cheek teeth. The Extraction Forceps can be disassembled and reassembled again in many differentways like the spreaders.

For horses with small heads we developed the "Pony Model". The actual version is dedicated to Professor Padraic Dixon from the University of Edinburgh who, at the end of the 20th century, made great efforts to reintroduce oral extraction techniques.



REVERSE FULCRUM

Thanks to this new Reverse Fulcrum Extension there is no extra Reverse Fulcrum Forceps needed. This system works both for lower and upper cheek teeth. It is especially recommended for 06 and 07's. It can be a quite helpful approach for lower 09 and 10 as well due to the curved nature of these teeth.

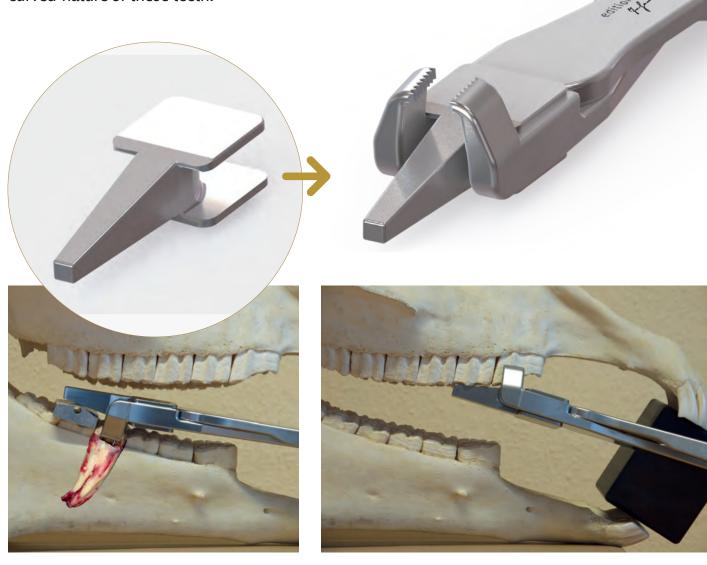

ArtNo.#	Description
PEGA_321020	Molar Forceps - Pony
PEGA_321001	Molar Forceps No. 1 – large
PEGA_321002	Molar Forceps No. 2 - large
PEGA_321010	Reverse Fulcrum Extension (fits Molar Forceps No. 1 and No. 2)

KIT FOR MINIMALLY INVASIVE TRANSBUCCAL EXTRACTION (MTE)

The MTE technique is indicated in the following cases:

- Force applicable to the crown is insufficient
- Crown is broken off or too fragile
- Malformation of the crown
- Removal of root fragments

To use the MTE Kit, it is necessary to attend a workshop where this technique is taught in-depth.



ArtNo.#	No.	Description
PEGA_322000		MTE Kit
PEGA_UM301140	1	Autoclavable container
PEGA_UM301090	2	Hammer
PEGA_UM301100	3	Dental chisel – plain
PEGA_UM301102	4	Dental chisel – offset
PEGA_UM312001	5	Short handle
PEGA_UM312010	6	Luxator
PEGA_UM312011	7	Sharp spoon
PEGA_UM312012	8	Loop
PEGA_UM301051	9	Drillingsleeve

ArtNo.#	No.	Description
PEGA_UM301039	10	Short blunt trocar
PEGA_UM301038	11	Short sharp trocar
PEGA_UM301014	12	Short trocar sleeve
PEGA_UM301041	13	Long blunt trocar
PEGA_UM301040	14	Long sharp trocar
PEGA_UM301015	15	Long trocar sleeve
PEGA_UM301130	16	Root fragment forceps
PEGA_UM301045	17	Extraction screw
PEGA_UM301044	18	Threading tap
PEGA_UM301042	19	Drill

DENTAL PICKS

DENTAL PICK SET - BASIC

This set contains a 8° handle, a set of 5 mm / 0.19 inches luxators in different lengths and a non-autoclavable screwdriver. It can be used to loosen gingiva, to remove tooth fragments or remaining roots. The luxators have a fine and sharp tip and can be applied in any turning direction.

The short and blunt insert is intended to fix the screw during transport.



ArtNo.#	Description
PEGA_310100	Dental Pick Set – Basic T–Handle– 8°, Wide Luxator Set, 5 mm (0.19 inches), Screwdriver
PEGA_310101	T-Handle - 0°
PEGA_310102	T-Handle - 8°
PEGA_310103	T-Handle - 18°
PEGA_310104	T-Handle - 30°
PEGA_310000	Wide Luxator Set -5 mm (0.19 inches)
PEGA_310010	Narrow Luxator Set -4 mm (0.16 inches)
PEGA_310020	Chisel Set
PEGA_310030	Probe Set
PEGA_310040	Sharp Spoon Set
PEGA BLOCK	Storage block, autoclavable (hard tissue)
PEGA_310106	Storage block, autoclavable (metal)

DENTAL PICK SYSTEM: 4 DIFFERENT HANDLES

Dental picks are intended for intraoral use.
There are 4 different handles, corresponding to 4 different angles:



The **0**° Handle can be helpful for the mandibular jaw.



The 8° Handle is the most often used one (up to 80% of the cases can be solved with it). Inserts can be used in a backward direction as well (helpful for lower premolars).



The 18° and 30° Handles are used for the more caudal teeth.

INSERTS

Additionally, there are different types of inserts available. They can all be used with the short straight handle as well to have even more versatile applications.



Narrow Luxator Set - 4mm/0.16 inches (5 lengths)
The narrow luxators are used for extra small fragments.



Chisel Set (4 lengths)
The chisels are used to get into the interdental space to loosen or remove tooth fragments.



Probe Set (5 lengths)
These probes are particularly useful for exploring deep fistulacanals and other dental and periodontal pockets.



Sharp Spoon Set (4 lengths)
These curettes are useful for cleaning the tooth socket and removing deep root fragments.

SHORT HANDLE EXTRACTION KIT




This kit consisting of a short handle with different luxators is very useful for the extraction of wolf teeth, incisors, canines and deeply located root fragments.

It contains straight, concave and convex curved luxators in a 4 and 5 mm / 0.16 and 0.19 inches version.

Particularly the concave curved tips allow a very good access to the palatal side of the wolf tooth and enable a careful extraction with minimal damage.


The short handle is designed to allow the gentle usage of a hammer as well. Furthermore, this short handle can be combined with all differentinserts from the Dental Pick System.

Shapes

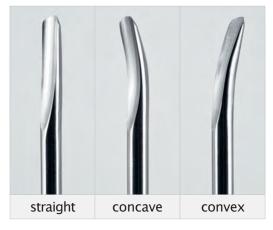

Widths



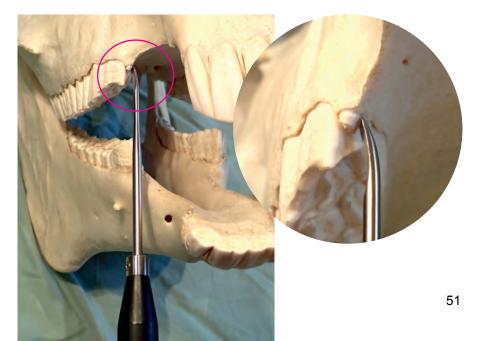

ArtNo.#	Description
PEGA_312001	Handle (short, straight)
PEGA 312004	Luxator Insert 150/4straight
PEGA_312002	Luxator Insert 150/5straight
PEGA_312005	Luxator Insert 150/4concave
PEGA_312003	Luxator Insert 150/5concave
PEGA_312006	Luxator Insert 150/4convex
PEGA_312007	Luxator Insert 150/5convex
PEGA_312040	Wolf Tooth Extraction Kit Contains a short handle and 4 luxators (straight and concave curved inserts in 4 and 5 mm width)
PEGA_312041	Extraction Kit with short Handle and 6 Luxators Contains a short handle and 6 luxators (straight, concave and convex curved inserts in 4 and 5 mm width)





HDE COMMERCIAL GUARANTEE

1. SCOPE

In addition to the Legal Guarantee, HORSE DENTAL EQUIPMENT (HDE) offerspractitioners an HDE Commercial Guarantee which applies to all new HDE-branded hand-held power equipment purchased by the practitioner from HDE or one of its official distributors*, from June 1st, 2016.

2. DEFINITION

a) Product Commercial Guarantee

HDE guarantees to the practitioner the breakage of his/her HDE-branded hand-held power equipment for a 12-month period following its date of purchase from HDE or one of its official distributors*. This guarantee applies to:

- -Broken parts (except "wear" parts -paragraph 3.b)
- -Labour costs
- -Transport costs.

b) After-SalesService Commercial Guarantee

All repairs made by HDE on HDE hand-held power equipment are guaranteed 6 months.

This guarantee covers:

- -Broken parts replaced in the last 6 months (except "wear" parts -paragraph 3.b)
- -Labour costs for the replacement of these parts
- -Transport costs.

c) Free annual check-up

After the initial 12-month period of guarantee, HDE or Official Technical Centers* offer one Free Annual Check-up to extend the life of your equipment. It includes the following:

- -Complete test of the equipment
- -Disassembly
- -Control of bearings, gears, seals, transmission cables, diamond disc, electrical connections, springs
- -Cleaning
- -Greasing
- -Reassembly

Labour costs are free of charge, transport costs are at the practitioner's expense.

This free check-up is valid once a year (calendar year).

Further to this check-up,in case repairs are required, spare parts, repair labour and transportation costs are at the practitioner's expense. In which case, an estimate will be sent for approval before any repair is done.

3. IMPLEMENTATION

a) Product registration

The product commercial guarantee, the after-salesservice commercial guarantee and the annual check-upwill only be implemented if HDE products are registered on HDE's website within 15 days following the purchase by the practitioner from HDE or one of its official distributors*.

b) Disclamer

The commercial guarantee applies to all HDE hand-held power equipment except:

A - Following wear parts:

- -Diamond coating (Flexxidisc, Polyfloat burrs, axle burrs...)
- -Gears (Flexxidisc and Polyfloat)
- B If the problem is caused by an accident
- C If the damage is caused by the user
- D If the userfails to comply with the instructions manual and/or HDE cleaning recommendations*
- E If the serial number of the product is removed or modified
- F If the product is not disassembledin an Official HDE Technical Center* or modified with non-HDE parts
- G If the product was not registered on HDE's website within 15 days after its purchase by the practitioner



SUMMARY

POWER EQUIPMENT —	page 4
PROFILING EQUIPMENT —————	page 7
FLEXXIDISC ————————————————————————————————————	page 14
POLYFLOAT	page 16
MANUAL RASPS ———————————————————————————————————	page 25
STOMATOLOGY ————	page 31
MOUTH-OPENERS-	page 46
ACCESSORIES —	
CLEANING RECOMMENDATIONS —————	page 56



#1 European manufacturer of equine dentistry instruments

SETS

SET LEVEL 1

SET_L1

- -N1 engine (complete unit)
- -Quickflex reinforced drive shaft with spring 1.20 m
- -Flexxidisc long curved twist
- -Stainless steel brush with holder
- -Instrument case













SET LEVEL 2

SET_L2

- -EVOLUTION engine (complete unit with jacket)
- -Quickflex reinforced drive shaft with spring 1.20 m
- -Flexxidisc long curved twist
- -Handpiece for the 06s
- -Stainless steel brush with holder
- -Trolley instrument case





Also available with WATERCOOLING

SET_L2IR











SET LEVEL 3

SET_L3

- -EVOLUTION engine (complete set) with jacket
- -Quickflex reinforced drive shaft with spring 1.20 m
- -Flexxidisc long curved twist
- -Polyfloat long straight twist (without burr)
- -Polyfloat applecore Ø 14
- -Handpiece for the 06s
- -Stainless steel brush with holder
- -Trolley instrument case





Also available with WATERCOOLING SET L3IR













SET LEVEL 4

SET_L4

- -EVOLUTION engine (complete set) with jacket
- -Quickflex reinforced drive shaft with spring 1.20 m
- -Flexxidisc long curved twist
- -Polyfloat long straight twist (without burr)
- -Polyfloat applecore Ø 14
- -Polyfloat occlusal burr
- -Polyfloat long versatile cylindrical burr
- -Polyfloat conical burr long model
- -Handpiece for the 06s
- -Handpiece for the 11s
- -Stainless steel brush with holder
- -Trolley instrument case





Also available with WATERCOOLING

SET_L4IR























EVOLUTION ENGINE

EVOLUTION ENGINE (completeunit with jacket)

MEV_HDEBVEU

Wireless unit: work wherever and whenever you want

Autonomy: powerful & long-lasting battery for a daily practice (up to 10 horses)

Safety: the electronic clutch protects the equipment

Comfort and versatility: optimum weight distribution on shoulders and waist



Delivered with all major international plugs (European, UK,

Australian and US)

EVOLUTION DRIVE SHAFTS

EVOLUTION QUICKFLEX DRIVE SHAFT

- -Reinforced and silent
- -Stainless steel quality
- -No shear pins



Length 1 m FLEX_R90SQE

Length 1.20 m FLEX_R120SQE

Length 1.60 m FLEX_R160SQE



A rapid fastening QUICKFLEX system adapted to all HDE burrs. Thanks to the integrated clutch in the control box, no need for shear pins.

EVOLUTION ACCESSORIES

SHEATH

Right: ACC_FD Left: ACC_FG



BATTERY CAR CHARGER MEV_CVBE

To charge the battery while travelling

TRANSFORMER

MEV_TEFR







PNEUMATIC EVOLUTION **FOOTSWITCH**

MEV_PPCE




This new model of pneumatic footswitch offers maximal safety. It can be used on a damp surface without being damaged. On/off function only.

WATERCOOLING

Watercooling is an option which can be adapted to the EVOLUTION engine. Conceived to reduce heating when working on overgrowths as well as minimizing dentine dust.

KIT WATERCOOLING KIT_KWCVMES

-Capacity: 1.25 L

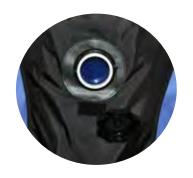



QUICKFLEX DRIVE SHAFT -WATERCOOLING

- -Quickflex reinforced drive shaft w/ spring 1 m for Evolution
- -Quickflex reinforced drive shaft w/ spring 1.20 m for Evolution
- -Quickflex reinforced drive shaft w/ spring 1.60 m for Evolution

FLEX_R90SQEIR FLEX R120SQEIR FLEX R160SQEIR



FLEXXIDISC -WATERCOOLING

-Flexxidisc long straight twist

-Flexxidisc long curved twist

FRA_GF7IR FRA_GF8TIR



POLYFLOAT -WATERCOOLING

-Polyfloat long straight twist (without burr)

-Polyfloat long curved (without burr)

FRA_GPOLYV2IR FRA_GPOLYCV2IR





N1 ENGINE

Powerful: for a regular equine dentistry practice

Simplicity: light and worn at the belt, it makes movements easy Adjustable speed: the practitioner adapts the speed to the treatment

Versatile: works with any type of HDE burr

N1 ENGINE (complete unit) MOT 1CREU

- -Power supply 220V or 110V
- -4000 -9000 rpm
- -Supplied with power cord + belt



European version (Swiss -UK versions available)

QUICKFLEX DRIVE SHAFT

- -QUICKFLEX system for a quick fastening, for all HDE drive shafts
- -Works with shear pins
- -Stainless steel
- -Available lengths: 1.20 m -1.60 m



1.20 m	
With stainless steel shear pins	FLEX_R120SQAI
With brass shear pins	FLEX R120SQAL



1.60 m		
With stainless steel shear pins	FLEX_R160SQAI	
With brass shear pins	FLEX_R160SQAL	

DRIVE SHAFT SHEAR PINS



STAINLESS STEEL SHEAR PINS
Bag of 20

Recommended with FLEXXIDISC & POLYFLOAT.



FLEX_LAIT

BRASS SHEAR PINS

Bag of 50

Recommended with all axle burrs.

N1 ENGINE ACCESSORIES

QUICK SHEAR PIN

FLEX_QA

- -Aluminum quality
- -The Quick Shear Pin allows to replace shear pins easily
- -Shear pins sold separately



CLUTCH

FRA_LCV3

- -The clutch is a protection system for the N1 engine
- -Weight: 0.165 kg



Thanks to the clutch it's no more necessary to use shear pins.



SHEATH

Right: ACC_FD Left: ACC_FG



PNEUMATIC EVOLUTION FOOTSWITCH

MOT_PPPC



N1 ENGINE FOOTSWITCH WITH SPEED CONTROL

MOT_006



63

HDE FLOAT

The HDE FLOAT is designed for occasional practice of equine dentistry.

HDE FLOAT-CURVED *

HDEF_220VCLC

-A full 360° rotation



HDE FLOAT-STRAIGHT*

HDEF_220VLC

-A full 360° rotation



FLEXXIDISC

The Flexxidisc is a versatile burr perfectly adapted for an efficient work on overgrowths as well as for routine profiling. It exists in different shapes and sizes adapted to every type of dentistry practice.

- -A high quality diamond disc (± 900 horses).
- -Diamond minimizes soft tissus damage.
- -The small height of the head eases work on the 11th.
- -The 360° rotating head allows optimal access (twist system on the curved version only).
- -The curved version follows the natural shape of the jaw, allowing better access in the back of the mouth.
- -Easy maintenance (no greasing necessary).
- -The ergonomic rubber handle allows an optimal grip.

FLEXXIDISC RANGE

Model	Curved	Straight	Length (cm)	Height head (mm)	Watercooling	Disc Ø (mm)	(1)
Long	~	~	54	27/12		32	
Standard	w orw/o extension	w orw/o extension	32	27/12		32	
Short		~	18	27/12		32	





FLEXXIDISC

FLEXXIDISC SHORT STRAIGHT

FRA_F7C

Length: 18 cm



FLEXXIDISC STANDARD STRAIGHT

FRA_F7

Length: 32 cm



FLEXXIDISC LONG STRAIGHT

FRA_GF7

Length: 54 cm



FLEXXIDISC STANDARD CURVED TWIST

FRA_F8T

Length: 32 cm



FLEXXIDISC LONG CURVED TWIST

FRA_GF8T

Length: 54 cm



The Polyfloat is a very versatile tool which allows a global dentistry treatment, from routine profiling to decay and diastema treatment, thanks to a complete range of 14 different burrs (the Polyfloat is delivered wihout burr).

It exists in different shapes and sizes, adapted to every type of dentistry practice.

- -The 360° rotating head allows optimal access (twist system on the curved version only).
- -The curved version follows the natural shape of the jaw, allowing better access in the back of the mouth.
- -A range of 14 different burrs for every type of treatment (profiling, decay, diastema, 06s).
- -Easy maintenance (no greasing necessary).
- -The ergonomic rubber handle allows an optimal grip.

POLYFLOAT RANGE

Model	Curved	Straight	Length (cm)	Watercooling
Long	~	~	54	
Standard	w orw/o extension	w orw/o extension	31	
Short		~	15	





<u>POLYFLOAT</u>

Non-contractual pictures: Polyfloat burrs sold separately.

POLYFLOAT SHORT STRAIGHT

FRA_POLYV2C Length: 15 cm



POLYFLOAT STANDARD STRAIGHT

FRA_POLYV2 Length: 32 cm



POLYFLOAT LONG STRAIGHT

FRA_GPOLYV2 Length: 54 cm



POLYFLOAT STANDARD CURVED TWIST

FRA_POLYCV2 Length: 32 cm



POLYFLOAT LONG CURVED TWIST

FRA_GPOLYCV2 Length: 54 cm



POLYFLOAT BURRS

The 14 HDE Polyfloat burrs benefit from a high quality diamond coating (± 900 horses). Diamond minimizes soft tissue damage.

The HDE range features 3 types of burrs:

APPLECORE BURRS

Applecore burrs are ideally used for a work on overgrowths, especially on the 11th (Ø14), 06s profiling (with the short straight Polyfloat) as well as final profiling. HDE offers a range of 6 burrs of different diameters and shapes.



DECAY TREATMENT BURRS

HDE offers 3 burrs of different sizes and shapes, specifically developed for decay treatments.



DIASTEMA TREATMENT BURRS

HDE offers 5 burrs of different sizes and shapes, for both lingual or vestibulary work.



11^s BURRS

Available in 3 different versions (straight, concave or convex), thoses burrs are specially designed to work on the 11^s. The compact, long and narrow head with a limited height (1 cm) facilitates an optimum accessin the back of the mouth.

- -The diamond coating allows to reduce easily the back hooks.
- -The concave and convex 11s burrs are adapted to the shape of the jaws for an easier access in the back of the mouth.
- -Easy maintenance (no greasing necessary).
- -Its rubber handle enables a perfect grip.



AXLE BURRS

MAXILLARY BURR

FRA_F1

-Length: 29 cm -Ø 15 mm / 5 mm

-Used with or without extension



-106 to 108 & 206 to 208:

Its cylinder shape allows fast & efficient profiling.

-109 to 111 & 209 to 211:

Its conical shape eases access to narrow areas between cheek and tooth.

Its conical protection applied to the cheek avoids contact with soft tissues and guides the burr on the tooth..

MANDIBULARY BURR

FRA_F2

-Length: 28 cm

−Ø 20 mm

-Used with or without extension



-306 to 308 & 406 to 408:

Its tapered sides ease access and 06° profiling.

-309 to 311 and 409 to 411:

Its protection has been designed to reach molars at the back of the mouth without damaging soft tissues.

ANGLE BURR

FRA_F3

-Length: 15 cm

-Ø 16 mm

-Used without extension



This burr is very versatile. Thanks to the tapered part of its diamond head, it can be used:

-For the 06s

-On incisors and hooks

BOWL BURR Ø 20

FRA_F4

-Length: 15 cm

-Ø 20 mm

-Used without extension



Its round shape offers impressive comfort and ease for the profiling of the 06s.

AXLE BURRS

BOWL BURR Ø 18

FRA_F4D18SP

-Length: 15 cm

- -Ø 18 mm
- -Without burr protection
- -Used without extension



DISC BURR

FRA_F5

-Length: 15 cm

-Ø 50 mm

-Burr thickness: 1 mm

-Protection thickness: 8 mm

-Used without extension



Used to reduce and realign incisor arcades.

6^S BURR

FRA_F6C

-Length: 16 cm

-Ø 18 mm

-Used without extension



Its curve allows an easy 06s shaping.

Perfect for the profiling of important hooks, thanks to its large diamond coating diameter.

ERGONOMIC HANDLE FOR 65 OR ANGLE BURR

FRA_PAFBS

- -Improves grip
- -Very easy to install
- -Fixing screw provided





This handle has been specially developed for practitioners looking for a better grip on their 06s (FRA_F6C) or angle burr (FRA_F3).

AXLE BURRS

VESTIBULARY BURR

FRA_F6

-Length: 29 cm -Ø 15 mm



Used with a slight opened mouth-opener it allows better comfort for the horse and the practitioner. Its less guarded protection allows a perfect working position.

NEW

PONY BURR

FRA_FP

-Length:43 cm

-Ø 9 mm



Specifically designed to treat ponies it has proven to be very useful on new world camelids as well. Its compact size and tiny head offer optimal accessin small mouths.

CANINE BURR

FRA_FSC

-Length: 14.5 cm -Weight: 152 g -Coil: Ø 6 mm

-Used without extension

.4.5 cm L52 g

Simple and easy profiling thanks to its miniature coil burr which is perfectly adapted for canines.

TOOTH SOCKET BURR

FRA_F9

-Length: 10 cm

-Used without extension



BOWL BURR Ø 8

FRA_F4D8

-Length: 16 cm

-Ø 8 mm

-Used without extension



ACCESSORIES

EXTENSIONS

Easy to setup, extensions make it possible to work without putting the hand in the horse's mouth while following the movement of the horse.

Recommended with the maxillary, mandibulary, vestibulary and medium Flexxidisc and Polyfloat burrs. 3 different versions are available.

SHORT EXTENSION

FRA_EXTC

-Length: 16 cm -Stainless steel

MEDIUM EXTENSION

FRA_EXTM

-Length: 22 cm -Stainless steel

LONG EXTENSION

FRA_EXTL

-Length: 28 cm -Stainless steel



STAINLESS BRUSH

-Fast and efficient cleaning of all diamond burrs

-Optional stainless steel holder

STAINLESS STEEL BRUSH

ACC_BNUE





ACC_BAS



INSTRUMENT CASE

HDE offers a range of instrument cases with a space dedicated for each piece of equipment.

TROLLEY INSTRUMENT CASE

RAN_MT

- -High shock resistance
- -Inner dimensions in cm (L x W x H): 52.1 x 28.7 x 19.1
- -Outer dimensions in cm (L x W x H): 55.9 x 35.6 x 22.9
- -Weight (empty): 5.2 kg



INSTRUMENT CASE

RAN_MALIV4

- -Dedicated space for each piece of equipment
- -Dimensions in cm (L x W x H): $53 \times 43 \times 22$
- -High shock resistant
- -Weight (empty): 4.4 kg



STANDARD INSTRUMENT CASE

RAN_MALR

- -Protective foam
- -Polypropylene
- -Dimensions in cm (L x W x H): 81 x 27 x 11
- -Weight (empty): 2.5 kg



POLYFLOAT BURR INSTRUMENT CASE

RAN_MFP

- -Protective foam
- -Polypropylene
- -Dimensions in cm (L x W x H): 23 x 19 x 6
- -Weight (empty): 0.26 kg





MANUAL RASPS







MANUAL RASPS

DIAMOND CAUDAL RASP


RAP_CAU

-Abrasive material: diamond D852

-Total length: 57 cm -Blade length: 50 mm -Blade width: 28 mm -Blade height: 3 mm

Rubber pistol handle for optimal grip.

Its narrow width allows easy access to the 11th caudal part.

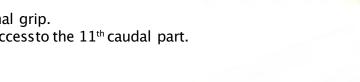

DIAMOND "S" RASP

RAP_SDIA

-Length: 17 cm-Width: 2 cm-Thickness: 0.5 cm

Essential for the finishing touch. With 2 different grains on one rasp.



TUNGSTENE CARBIDE "S" RASP

RAP_STUN

-Length: 30 cm -Width: 2.5 cm



3 HEAD RASP

RAP_TRI

-2 handles: 30 cm and 48 cm-3 tungsten heads: opened, closed

and straight angles



PONY FLOATHANDLE

RAP_BP1

- -Premolar pony float, 20° angle -38 cm
- -50 mm blade sold separately
- -20 cm neoprene handle
- -Stainless steel

Especially designed for ponies with narrow mouths.



77

MANUALRASPS

A range of quality rasps which minimize softtissuedamages and make interventions easy.

FLOAT HANDLE (without blade)

- -20 cm smooth rubber handle
- -Easy to clean
- -Titanium or tungsten carbide blades sold separtely -Stainless steel

-Stainless steel	Screw-on blades	Stick-on blades
Straight incisor rasp -30 cm 80 mm blade	B10	
Premolar rasp, 30° angle -38 cm 50 mm blade	B11	
Premolar rasp, 20° angle -43 cm 50 mm blade	B12	B22
Premolar or molar rasp, 10° angle -53 cm 80 mm blade	B13	B23
Straight lingual & buccal rasp -58 cm 80 mm blade	B14	B24
Straight lingual & buccal rasp -63 cm 80 mm blade		B25
Upper back molar rasp (upturned head), 15° angle -63 cm 50 mm blade	B16	B26
6s premolar roller rasp -43 cm 50 mm blade	B17	
Curved shaft 6s premolar roller rasp -43 cm 50 mm blade	B18	

BLADES & INSERTS FOR RASP HANDLES

TUNGSTEN BLADES

-Tungsten carbide blade 50 mm (to screw) REC_V5TU -Tungsten carbide blade 50 mm (to stick-on) REC_C5TU

-Tungsten carbide blade 80 mm (to screw) REC_V8TU

-Tungsten carbide blade 80 mm (to stick-on) REC_C8TU



TITANIUM BLADES

-Titanium	blade	50 mm	(to screw)	REC_V5TI
-Titanium	blade	50 mm	(to stick-on)	REC_C5TI
-Tungsten	blade	80 mm	(to screw)	REC_V8TI
-Tungsten	blade	80 mm	(to stick-on)	REC C8TI



TRIANGULAR INSERTS

-Tungsten insert for triangular rasp 3 grooves (unit)
-Tungsten insert for triangular rasp 4 grooves (unit)
-Titanium insert for triangular rasp 3 grooves (unit)
-Titanium insert for triangular rasp 4 grooves (unit)
REC_T3TU
REC_T4TI



ACCESSORIES

BAG FOR LONG RASP HANDLES

-Bag for long straight rasp 60 cm RAP_PRDL -Bag for short straight rasp 56 cm RAP_PRDC



STANDARD INSTRUMENT CASE

RAN_MALR

-Length: 81 cm -Width: 27 cm -Height: 11 cm -Protective foam

Useful to carry and protect rasps.



MANUAL RASP HOLDER

ACC_SUPPR





STOMATOLOGY





BURGESS KIT WITH OPEN INSERTS

BUR_KIT

- -4 open inserts (Ø 6 mm, 8 mm, 10 mm, 12 mm)
- -1 handle
- -1 curved extension
- -1 straight extension
- -1 rod (inserts)
- -1 rod (extension)
- -1 Burgess instrument box



BURGESS KIT WITH CLOSED INSERTS

BUR_KITF

- -4 closed inserts (Ø 6 mm, 8 mm, 10 mm, 12 mm)
- -1 handle
- -1 curved extension
- -1 straight extension
- -1 rod (inserts)
- -1 rod (extension)
- -1 Burgess instrument box



STOMATOLOGY

KIT BURGESS (to compose)

BUR_KITM

- -1 handle
- -1 curved extension
- -1 straight extension
- -1 rod (inserts)
- -1 rod (extension)
- -1 Burgess instrument box



Compose your own kit by choosing from the following inserts:

CLOSED INSERTS

6 diameters (4 mm, 6 mm, 8 mm, 10 mm, 12 mm, 14 mm)



OPEN INSERTS

4 diameters (6 mm, 8 mm, 10 mm, 12 mm)



31



KIT LECLAIR

KIT_EXT

- -Wolf tooth elevator n°1
- -Wolf tooth elevator n°2
- -Straight elevator
- -Bone curette 8 mm
- -Gouge forceps 45° with spring
- -Curved speculum 17 cm
- -Curved speculum 45° open 20 cm
- -Instrument box $(25 \times 15 \times 5 \text{ cm})$



Everything you need for an easy and risk-freewolf tooth extraction.

WOLF TOOTH ELEVATORS

- -5 models: curved right, curved left and straight
- -Stainless steel







WOLF TOOTH ELEVATOR

3 vestibular elevators

3 lingual elevators



STO_EDL20V3.5



STO_EDL20L3.5



STO_EDL20V5.5



STO_EDL20L5.5



STO_EDL20V7.5



STO_EDL20L7.5

Also available as a set of 6: STO_S6EDL





SMALL FORCEPS

-Stainless steel

SMALL FORCEPS

Straight forceps closed



Straight forceps open



DAV_DBO

90° forceps closed



DAV_C90BF

90° forceps open



DAV_C90BO

Curved forceps 17 cm



DAV_CTC

45° FORCEPS

DAV_DBF

DAV_C45BO

-Length: 20 cm -Stainless steel

Its 45° head makes the wolf tooth extraction easy without soft tissue damage.

Its length guarantees optimal comfort during extraction.



90° FORCEPS

-Length: 25 cm -Stainless steel



SPECIAL LEFT DAV_2590SPEWG

SPECIAL RIGHT DAV_2590SPEW





GOUGES & CURETTES

BONE CURETTE

-Stainless steel

-Bone curette 6 mm STO_CUR6
-Bone curette 8 mm STO_CUR8

-Bone curette 10 mm STO_CUR10







GOUGE FORCEPS STO_PIN



Its 45° angle makes access to the soft tissue around the wolf tooth easier. Allows not to be bothered by the lips nor by the mouth opener straps.

35



GOUGES & CURETTES

BONE RONGEUR

STO_GOU

-Length: 29 cm

-Head diameter: 10 mm



GOUGES

- -Stainless steel
- -2 types (straight and curved)
- -Diameters: 6 8 10 mm

GOUGES STRAIGHT

GOUGES CURVED





PICS & SONDS


DENTAL PIC

-Length: 46.5 cm -Height: 4 cm

-Width of the stainless steel handle: 8 cm



Its better grip gives more strength when rotating. The handle acts like a visual and tactile cue while working in the horse's mouth.

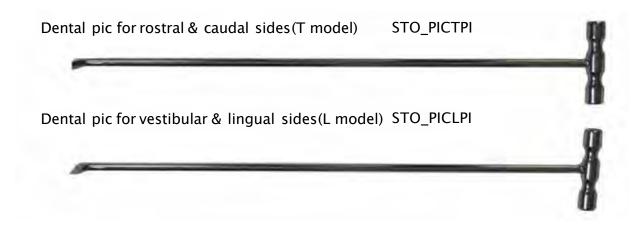

DENTAL SOND

-Length: 45 cm

-2 models: left and right

-Ideal for preparing molar extractions

-Stainless steel



STO_SONG



PICS & SONDS

EXPLORER SOND

STO_SE

- -Allows to probe and diagnose periodontal tissues
- -Total length: 41 cm
- -Stainless steel



DENTAL HOOKS

CLASSIC DENTAL HOOK

STO_CD

- -Very resistant
- -Easy to disinfect
- -Length: 47 cm
- -Height of the hook: 2 cm
- -Stainless steel

DENTAL HOOK -SMALL MODEL

STO_CDF

- -Total length: 41 cm
- -Diameter of the hook: 2 cm
- -Stainless steel





Allows to curette vestibular parts before extraction and to eliminate interdental residues.



LARGE FORCEPS



Its ergonomic handle and adjusted head ease the extraction of milk teeth.



This instrument loosensthe molar to extract by leaning on the adjacent molars.



MOLAR FORCEPS

MOLAR FORCEPS WITH SCREW

DAV_MV

-Total length: 53.5 cm

-Possible height: from 1.5 to 3.5 cm

-Stainless steel



Thanks to its adjustable screw which acts like a lever, this forceps ensures better comfort during molar extraction.

MOLAR FORCEPS WITHOUT SCREW DAV_M

-Total length: 53.5 cm

-Stainless steel



FULCRUM

-Length: 43 cm

-Stainless steel handle brass inserts

-2 different models



Fulcrums are used as a lever for the molar forceps. During extraction, the different sizes and shapes will help remove the molar without difficulty.



FRAGMENTS FORCEPS

4-ROOTMOLAR FORCEPS

DAV_F4ML

-Length: 53.5 cm -Weight: 1.8 kg -Stainless steel

Easily extract molars thanks to its specially-designed angles.



3-ROOTMOLAR FORCEPS

-Length: 53.5 cm -Weight: 1.8 kg

-2 different models (right and left)

-Stainless steel



RIGHT MODEL

DAV_F3B



LEFT MODEL

DAV F3BG



Its three-pointgrip is perfect for removing triangular fragments.

The side with the two points is positioned on the left or the right side of the tooth (depending on the forceps) and the other on the fractured side.

2-ROOTFRAGMENT FORCEPS

DAV_F2M

-Length: 53.5 cm-Weight: 1.8 kg-Stainless steel



A forceps with 2 long teeth for small molars fragments.

The diameter and length of its teeth capture and remove fragments within the tooth socket.

DENTAL MIRROR WITH POUCH

STO_MI

-Length: 47.5 cm -Diameter: 50 mm -Angle: 25°

-Angle: 25° -Anti-mist mirror -Stainless steel



DENTAL MIRROR -SMALL MODEL

STO_MDD40

-Length: 41 cm -Diameter: 44 mm

–Angle: 25°

-Anti-mist mirror -Stainless steel



ACC_BI

-Dimensions in mm (L x W x H): $250 \times 150 \times 50$

-Weight (empty): 0.7 kg

-Stainless steel



WOLF TOOTH NERVE BLOCK SYRINGE -25 cm *

-Syringe for cartridges 1.8 ml ACC_SAL251.8 -Syringe for cartridges 2.2 ml ACC_SAL25



Ideal for administering local anesthetic for nerve block prior to wolf tooth extraction. (Needle and local anesthetic not included)

GENERAL / MANDIBULAR NERVE BLOCK SYRINGE -40 cm *

-Syringe for cartridges 1.8 ml

ACC_SAL401.8

-Syringe for cartridges 2.2 ml

ACC_SAL40



Ideal for administering mandibular local anesthetic for nerve blocking and general nerve block prior to extraction. (Needle and local anesthetic not included)

MALLET

STO_M

-Length: 16 cm-Diameter: 4 cm-Weight: 0.9 kg



94

SET OF 3 TREPHINES

STO_K3T

-3 diameters: 19 -22 -25 mm

-Stainless steel



MOUTH-OPENERS-ACCESSORIES







MOUTH-OPENERS

OPTIMUM MOUTH-OPENER (WITHOUT PLATES)

OB_OPT

-Weight: 3.4 kg (with plates)

-Opening: 15.5 cm -Biothane straps -Stainless steel



Safety & comfort: 29 safety catches and wide arms ensure optimum comfort for the patient and practitioner.



MOUTH-OPENERS

STANDARD MOUTH-OPENER (WITHOUT PLATES)

OB_STA

-6 safety catches -Weight: 3 kg

-Biothane straps

-Stainless steel

Also suitable for young horses up to 2 years old.



PONY MOUTH-OPENER (WITHOUT PLATES)

OB_PON

-5 safety catches

-Weight: 2.3 kg

-Opening: 7 cm

-Biothane straps

-Stainless steel



PONY MOUTH-OPENER PLATES

OB_PLAP

(Sold separately)



INCISOR MOUTH-OPENER OB_INC

-Flexible PVC pipe -Adjustable strap



OPTIMUM & MILLENIUM ARCH SPREADER OB_ARCE



- -For incisor work
- -Easy to install and very confortable for the horse

LATERAL MOUTH-OPENER

OB_LAT

-Weight: 0.4 kg -Length: 30 cm

-Diameter of the nylon piece (small): 3.5 cm -Diameter of the nylon piece (large): 4.5 cm



CHEEK RETRACTOR

OB_RJ

-Weight: 0.25 kg

-Height: 17 x 10 x 5 cm



For better visibility during diagnosis before and after the intervention.

47

98



Our plates can be used with Optimum, Millenium & Standard mouth-openers-Stainless steel.

CLASSIC PLATES

Standard OB_PLAC







Like the arch spreader, these plates make work on the vestibular part easier.

INCISOR MALOCCLUSION PLATES

Standard OB_PLABC

Large OB_PLABL





Large surface tooth plates for an easy access on horses with incisor malocclusions.

INCISOR PLATES OB_PLAI



Specially designed for safe work on the incisors.

SURGICAL LAMP HEADLAMP WITH CARRYING CASE

ACC_LCV3

- -Autonomy up to 4 hours
- -Rechargeable battery (charging time around 4 hours)
- -Working distance: 420 mm
- -LED: 3 W / 60 000 lx
- -Luminous intensity: 5500 Kelvin
- -LED lifetime (50 000 hours)
- -Lamp's weight: 280 g
- -Delivered in a carrying case



This surgical headlamp is easy to use, lightweight and wireless. Adjustable height and focusable lamp which ensures optimal lighting in the horse's mouth.

MAGNETIC LAMP

OB_LOBV5

- -Integrated magnetic holder
- -Enhanced connections
- -Battery autonomy: 12 hours
- -Light charge: 130 lx
- -Angle: 170°



RECHARGEABLE LED HEADLAMP

ACC_LFLV2

- -1 washable elastic headband
- -1 battery, 1 power adapter, 1 USB charger
- -Integrated battery status by color system
- -Autonomy up to 60 hours
- -Charging time: 5 hours
- -Luminous intensity up to 300 lx
- -Total weight: 165 g



The new LED rechargeable headlamp is brighter, smarter, more durable and comfortable than the "standard" H7R version.

Can be used with 4 alkaline batteries AAA.

EQUINE DENTALSCOPE

NEW!

Equine dentalscopes are becoming increasingly useful in equine dentistry. They enable a more thorough examination of the oral cavity & dental arcades and allow a better detection of dental pathologies and utilise the camera and video capabilities of your laptop, phone or tablet.

BENEFITS

- Clear and precise diagnosis
- > Better patient treatment and follow-up
- > Better client dental education
- Added value and service









Starter Set-1

-Dentalscope for connection to a display device with cable

ACC_CD

Starter set + Wireless USB Connector-2

Dentalscope with USB wireless connector for phone / computer /tablet.

ACC CDUSB

Starter Set + Wireless Monitor-3

Dentalscope with a 10'' wireless monitor including recording system, micro-SD card and storage case.

ACC_CDM

FEATURES

- > 3 monitoring options: TV screen (option 1), Wireless USB connector for phone / tablet / computer (option 2), Wireless screen (option 3)
- Available on PC Win 7, Win 10, Android phones and tablets
- ➤ 692 mm total length
- ► Head dimensions: width 17mm x thickness 17,2mm
- Rounded head on top
- > Stainless steel
- Weight: 360 g

50

RINSING PEAR

ACC_POIR

-Capacity: 250 ml

RINSING PISTOL

ACC_PR200ML

-Capacity: 200 ml





RINSING SYRINGE

ACC_SERIN

-Capacity: 400 ml

ACC_SR600ML

-Capacity: 600 ml



STAINLESS STEEL BUCKET

ACC_SI12L

-Capacity: 12 L

-Protective rubber layer on the bottom of the bucket

-Stainless steel



DENTAL CHART

ACC_BSDF

- -100 dental sheets (50 sheets of carbon paper)
- -All necessary dental care on a detailed chart

Also availalble:

In Spanish ACC_BSDESP
In French ACC_BSDF
In German ACC_BSDALL



«LA DENTISTERIE EQUINE» BOOK

ACC_LDE

-Authors: Holger Schlemper -Isabell Herold

Tilman Simon
-Size: 21 x 30 cm
-Pages: 184

-Publication year: 2014 -Language: French

-Editor: "Point Vétérinaire -Collection Atlas"



DENTAL PASTE

STO_PDD

- -Dental paste for diastema filling
- -Includes 2 components to mix

MAGNETIC HOOK ACC_CA

-Stainless steel -Max: 10 kg





WATER HEATER

ACC_THER

Very useful in wintertime to keep hands warm and fingertips sensitive.



TUNIC

-Size 1 to 6: ACC_TT1 to ACC_TT6

-Customized embroidery upon request



OSCILLATING HEAD STAND

ACC_STO

-Weight: 5.45 Kg

-Height min-max: 111-180cm -Width of the headrest: 50 cm

-Aluminum

-Washable and removable cover





Thanks to its oscillating system, the column follows the horse's moves, avoiding any risk of injury.

53

104

CLEANING RECOMMENDATIONS

You have purchased quality professional equipment. Please find hereafter some advice on how to optimize its use and longevity.

GENERAL RECOMMENDATIONS

- -Make sure each aspect of your equipment is always clean & dry before storage.
- -HDE recommends at least 1 annual overhaul (or every 500 horses) of your complete equipment in an approved technical center*.
- -HDE does not cover damages caused by the use of bleach.
- -Find all of our recommendations in our products'instruction manuals.



DIAMOND BURRS

- 1 -Start the engine with the handpiece attached.
- 2 -Run it at high speedwith the head pointing down to remove any remaining water.
- 3 -Clean the diamond burr by running it over the HDE stainless steel brush.
- 4 Continue until the remains of dentine have completely disappeared.
- 5 -Clean with a cloth moistened with clean water and a disinfectant.
- 6 -Thoroughly dry with a clean, soft cloth.
- 7 -Store the burr in its dedicated place within your clean, dry case.

The burr is resistant to water projection, do not fully immerse in water. Always use clean water without corrosive nor abrasive products. Only use the HDE stainless steel brush.

DRIVE SHAFT

- 1 -Clean the outside sheath with a damp sponge.
- 2 -Dry thoroughly with a clean, soft cloth.
- 3 -When finished working, detach the drive shaft from the engine and the handpiece.
- 4 -Store the drive shaft in your clean, dry case.





Ensure not to bend the drive shaft in an excessive radius of curvature.

ENGINE

- 1 -Turn off the engine before cleaning.
- 2 -Clean the engine with a damp sponge.
- 3 -Thoroughly dry the engine with a clean, soft cloth.
- 4 -Store the engine in your clean, dry case.

HDE recommends to use the engine at medium speed (speed 3).

CASE

- 1 -Clean the case with a damp sponge.
- 2 -Thoroughly dry the case with a clean, soft cloth.
- 3 -Store in a clean and well ventilated place.

Make sure the case is opened from time to time, to ensure your equipment is well ventilated.



MAI Animal Health

NZ Equine Dental Catalogue



Treating Animals. Differently

AlumSpecTM Poental Products uc PATENTED

- Constructed of lightweight aluminum, netweight of 6 lbs
- Durability is not compromised by light weight design -able to with stand bending forces in dental procedures
- Use of replaceable synthetic straps in the "cinch-lock" locking systemallows for infinite degrees of access to the patient's oral cavity
- Easy to dean and is pliable at all temperatures
- May be used on all sizes of horses from miniature to draft
- Recommended for use with sedated animals for the safety of the handler and horse

ltem	Product Description	UOM
53600	AlumiSpec™ Equine Speculum	Each
53602	Replacement OpeningLatch Strap for AlumiSpec™ & X Spec™	Each
53603	CamBuckle for Opening Latch Strap for AlumiSpec™ & X Spec™	Each
53676	ReplacementPoll Strap for AlumiSpec™&X Spec™	Each
53675	ReplacementBolt Kit for AlumiSpec™&X Spec™, 4 Bolts	Each



X Spec[™] Equine Speculum, Wide Poental Products.

PATENTED

- Next generation equine speculum, an enhangement of the original AlumiSpec™, known and trusted throughout the industry
- Wider design provides morespace between the horse's cheek and speculum to accommodatedicital radiography plates for obtaining dental films in a variety of positions
- Constructed of aircraft grade aluminum, CNC milled precision pieces, TIG welded with a durable powdercoat finish
- Use of replaceable synthetic straps in the "cinch-lock" locking systemallows for infinite degrees of access to the patient's oral cavity
- Maybe used on all sizes of horses from miniature to draft
- Recommended for use with sedated animals for the safety of the handler and the horse
- Tested on hundredsof patients in an equine dental practice
- Netweight6lbs

ltem	Product Description		UOM
53605	X Spec™ Equine Speaulum		Each
53602	Replacement OpeningLatch Strap for AlumiSpec™ & X Spec	гм	Each
53603	CamBudde for Opening Latch Strap for AlumiSpec™&XSp	ec™	Each
53676	ReplacementPoll Strap for AlumiSpec™&X Spec™	ltem	Each
53675	ReplacementBolt Kit for AlumiSpec™&X Spec™, 4 Bolts	53425	Each



VETERINARY Dental Products to VETERINARY

MagneticSpec BiteLite™

VETERINARY Dental Products

- Fits moststyles and sizes of speculums by attaching to the upper bite plate
- Contains 2 LED lights
- Uses 3"AA" replaceable batteries
- Manufactured of heavy, durable plastic and is water resistant
- Metal protective plate protects from possible equipment damage
- Magnetic light uses 4 neodymium magnets with a combined rated 26.49 lbs of pull force, two 1/2" at 1663 lbs and two 3/8" at 9.86 lbs
- Coiled cord adds moreflexibility
- Mountingofspeculumplate is available upon request

Item	Product Description	UOM
53425	MagneticSpec BiteLite™	Each
53425-2	Replacement Battery Box	Each

Magnetic Speculum Plate Poental Products

- Magnetic plate designed for use with non-magneticspeculums
- Enables use of Magnetic Spec BiteLite™ with non-magnetic speculums

1		D 1 (D 12	
	Item	Product Description	UOM
	53425-1	Magnetic Speculum Plate	Each



OtoKleen



OtoKleen™Ear Cleansing Solution for Dogs&Cats (Moi)

Otokleen™ is designed as a non-irritating surface deaner to lift ear canal dirt and debris from the skin surface and removed via the vacuum action of the Auri-Kleen™ Irrigation / Aspiration Machine.

- Powerful waxdissolver
- Non-initating formula
- Can be used with the Auri-Kleen™lrrigation /Aspiration Machine
- Clean several ears with one bottle

Item	Product Description	UOM
48103	OtoKleen™ Ear Cleansing Solution for Dogs & Cats -16 oz	Each



Ahighquality affordable set of Surgical Loupes with 2.5X magnification and 19" working distance. Surgical Loupes are ideal fordental and surgical procedures. Adjustable inter-ocular barrels allow for multiple users.

Kit indudes:

- Loupes
- Head strap
- Safety frames
- Tool kit for adjustments
- Carryingcase

Item	Product Description	UOM
48320	Surgical Loupes, 2.5X	Each



LED Clip Light for Surgical Loupes WMDS-VET

- 100 lumens of bright white light
- 100,000hourLED life expedancy
- Operates on 3 "AAA" batteries (included)

Item	Product Description	UOM
48315	LED Clip Light for Surgical Loupes	Each
48314	Power Cord, 10' with 110-VoltTransformer&Switch for #48315Clip Light	Each



5-WattLED Headight WMDS-VET

- Extra-brightwhite, 5-WattLED
- 300 lumens of LED light for 100,000 working hours
- Fully adjustable headbandand overhead strap
- Adjustable light spot size, brightness and position
- 2 rechargeable batteries with charger
- Optional: Compatible with 2.5x Surgical Loupes


Item	Product Description	UOM
48312	5-WattLED Headlight with 2 Rechargeable Batteries	Each
48322	2.5x Magnification Surgical Loupes for 5-WattLED Headlight	Each









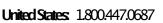





Dual LED Headight Poducts...

- Features two slim-lineLED spotlights mounted vertically following the nose line
- Lightweight headlight offers exceptional brightness
- Fully adjustable headbandcomfortablyfits all head sizes
- Mountedbattery boxrequires 3"AA" batteries (included)
- Battery boxis removable for convenience of changing batteries or replacement.

ltem	Product Description	UOM
53422	Dual LED Headlight	Each
53425-3	Replacement Battery Box for Dual LED Headlight	Each



LED Headlight



- Long-lastingLED bulb has three brightness settings
- Fully adjustable headbandcomfortablyfits all head sizes
- Brightness settings
 - ♦ High =1,000 Lumens, approximately 2 hour runtime
 - ♦ Medium=600 Lumens, approximately 3 hour runtime

Indudes:

- 2 removable, rechargeable battery packs
- Battery extensioncord
- Wall charger

Item	Product Description	UOM
53440	LED Headight	Each



Single LED Headlights WMDS-VET

The Single LED Headlight is a battery-powered LED headlight with 100,000 hour LED (light emitting diode) for surgical procedures or general examination. Two-wayheadbandadjustment on this light weight headlight makes it the ideal portable lighting solution.

- Overhead strap for comfort
- Adjustable headband
- Adjustable light position
- 100 lumens of bright white light
- Operates on 3 "AAA" batteries (#48305) (included) or "AA" batteries and 110-Voltransformer (#48316) (included)
- Optional, interchangeablelenses available

ltem	Product Description	UOM
48305	Single LED Headlight	Each
48316	Single LED Headlightwith "AA" Battery Pack & Transformer	Each



Barrhart Curette 5/6 USI BSI

	and Institute	
Item	Product Description	UOM
54550	Barnhart Curette 5/6	Each



McCals Curette 17/18 USI BSI

Item	Product Description	UOM
54551	McCalls Curette 17/18	Each



Columbia Curette 13/14 USI BSI

Item	Product Description	UOM
54552	Columbia Curette 13/14	Each



Gracey Curette 13/14 USI BSI

Item	Product Description	UOM
54557	Gracey Curette 13/14	Each



Acid Etch Cel (NOVADENT)

Item	Product Description	UOM
54604	Acid Bitch Cel, 12 ml	4/Pack



ROLINK* used Devel Afficial grant are developed afficial gran

ProFi™ Howable Composite

ltem	Product Description	UOM
54605	ProFil™ Flowable Composite,1 gram	4/Padk

ProLink™ Bonding Agent

Ite m	Product Description	UOM
54606	ProLink™ Bonding Agent, 5 ml	Each

Mixing Spatula (NOVADENT)

Item	Product Description	UOM
54609	MixingSpatula	100/Pack

Micro Applicator Brush (NOVADEDT)

	• •	
Item	Product Description	UOM
54610	Micro Applicator Brush -Super Fine	100/Padk

Flexi-Float Ultra Complete Kit Dental Products PATENTED

The Flexi-Float™ Ultra is a lightweight, compact, all-in-one, battery-powered floatsystem. The slender design and movable disc angle allows practitioners to maneuver around the dental arcade, along with its shaft which rotates 180 degrees for use on both the upper and lower arcade. This well balanced system has fully integrated irrigation, suction and trigger speed control which allows for quick removal of hooks in the mouth. The Flexi-Float" Ultra is battery powered making it completely mobile. The rechargeable battery lasts between 8-10 horses. The diamond discs are smoother and last longer than carbide discs and can be quidkly changed out in the field, taking less than one minute. Clean daily and grease after every 100 uses. Aluminum carrying tote sold separately.

Kit indudes:

- Flexi-Float™Ultra
- Battery-Powered Vacuum
- Li-Hon Batteries for Vacuum(2)
- Battery Chargerfor Vacuum
- **ImigationSystem**
- MaintenanceKit (Grease Gun, Brush)
- 14.4-VoltNiMHRechargeable Batteries (2)
- Battery Charger

Protective Carrying CaseInstruction Manual		
Item	Product Descritpion	UOM
53050R	Flexi-Float™ Ultra CompleteKit with RemoteWired Battery	Each
53053	RemoteWired Battery for Flexi-Float™ Ultra & Versa-Float™	Each
53054	Battery Charger for Flexi-Float™ Ultra & Versa-Float™	Each

Hexi-Float™ Ultra Aluminum Carrying Tote Dental Products...

Optional tote to carry Flexi-Float™ Ulfra instrumentandaccessories

	Product Description	UOM
53051	Flexi-Float™ Ultra. AluminumCarrying Tote	Each

Disc, 1" Diameter, 1/4'Shank, Diamond Dental Products...

- Diamondchip reduces tissue clamage
- 1"diameter with 2 holes for easy replacement
- 1/4"shank
- Fits Flexi-Float™ Ultra and Versa-Float™ Large Head Shaft

ltem	Product Description	UOM
53037	Disc, 1"Diameter, 1/4"Shank, Diamond	Each

53037

Diasterna Burs Poental Products

Fits Flexi-Float™ Ultra and Versa-Float™ Large Head Shaft

Item	Product Description	UOM
53038	Diastema Bur #1, Diamond	Each
53039	Diastema Bur #2, Diamond	Each
53040	Diastema. Bur #3, Diamond	Each









Versa-Float** Kits **P Dental Products... PATENTED

The Versa-Float" Kit is the latest advancementin Equine Dentistry, providing practitioners a mobile instrument to enhance floating techniques. The Versa-Float" dental systemoleivers a newstandard to power dental instruments. The variable speed unit is controlled by the trigger on the pistol-grip handle, allowing the practitioner more control and accuracy in the mouth. Kit customization available.

Kit includes:

- Versa-Float™ PlatformMotor
- Versa-Float™ Large Head Shaft with 1"DiamondDisc
- Versa-Float™ Small Head Shaft with 5/8"DiamondDisc
- Irrigation Systemfor the Large Head Shaft
- Irrigation Systemfor the Small Head Shaft
- MaintenanceKit
- 14.4-VoltNiMHRechargeable Batteries (2)
- Battery Charger
- Protective Carrying Case
- Instruction Manual



Item	Product Description	UOM
53079R	Versa-Float™ Kit, Large &Small Head Shafts with Platform& Irrigation with RemoteWired Battery	Each
53075R	Versa-Float" Kit, Large Head Shaft with Platformwith RemoteWired Battery	Each
53076R	Versa-Float™ Kit, Large Head Shaft with Platform& Irrigation with RemoteWired Battery	Each
53077R	Versa-Float™ Kit, Small Head Shaft with Platformwith RemoteWired Battery	Each
53078R	Versa-Float™ Kit, Small Head Shaft with Platform& Irrigation with RemoteWired Battery	Each
53028	Versa-Float™ Large Head Shaft	Each
53029	Versa-Float" Irrigation Attachmentfor Large Head Shaft	Each
53030	Versa-Float™ Small Head Shaft	Each
53031	Versa-Floati" Irrigation AttachmentforSmall HeadShaft	Each
53053	RemoteWired Battery for Flexi-Floati [™] Ultra & Versa-Floati [™]	Each
53054	Battery Charger for Flexi-Float™ Ultra & Versa-Float™	Each

Burs & Discs for Versa-Float™ Dental System-Small Head Shaft ✓ Dental Products ...

ltem	Product Description	UOM
53012	Diastema.Bur#1for Versa-Float™ Small, Diamond	Each
53013	Diastema Bur #2 for Versa-Float™ Small, Diamond	Each
53014	Diastema. Bur #3 for Versa-Floati™ Small, Diamond	Each
53015	Diastema Bur#4for Versa-Float™ Small, Diamond	Each
53020	Diamond Disc, 5/8'	Each
53021	Rotozip Bur, 1/8'Shank, 7/8'LOC, 2"OAL	4/Padk
Versa-Fl	oat"Small HeadShaft Accessories	
53032	Versa-Float™ Replacement Head, Small	Each
53033	Versa-Float™ Replacement Collet forSmall Head	Each
53044	Bur#2,Style DSpecial, FG, 1/8'Shank, 1/2'LOC, 3"OAL - Diamond	Each
53044-1	Bur#2,Style DSpecial, FG, 1/8'Shank, 1/2'LOC, 3"OAL-Diamond	2/Padk



ltem	Product Description	UOM
53037	Disc, 1"Diameter, 1/4"Shank, Diamond	Each
53038	Diastema Bur #1, Diamond	Each
53039	Diastema Bur #2, Diamond	Each
53040	Diastema Bur #3, Diamond	Each










De-MaX^M Reciprocating Float Dental Products

The De-MaX"III Reciprocating Float is a reliable reciprocating systemutilizing DeWalt® technology.

- Three different handle angles for ease in use
- Reciprocating open has been modified from a 9/16'stroke to a 3/8'stroke
- Low profile allows the user complete maneuverability around the dental arcade
- Excellent for removing caudal hooks and for finishing upper and lower molars
- 12-VoltLi-ion battery
- Kit indudes:
 - 19' Float Shafts (2); 1 Straight and 1 Angled
 - Batteries (2)
 - Charger
 - Fine Radius End LowProfile Blades (2)

Item	Product Description	UOM
53002	De-MaX*III Reciprocating Float	Each
53003	De-MaX"III Straight Shaft Float, 19"WorkingLength	Each
53003S	De-MaX"III Straight Shaft Float, 13"WorkingLength	Each
53004	De-MaX™II Angled Shaft Float, 19"Working Length	Each
53004S	De-MaX"III AngledShaft Float, 13"WorkingLength	Each

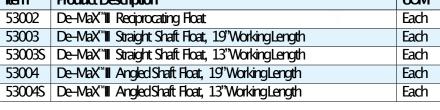

Power Handle™ Rotary Floats & Accessories ∨ VETERINARY Dental Products.

The Power Handle™ instrumentis an economical burpowersystemavailable in three sizes for accessibility to the incisors, premolars and molars.

- Cylindrical burmodelsindude a removable 2'soft tissue quard which protects the horse's cheek
- Constructed of heavy duty stainless steel
- 5/8'x2" carbide bur, solid cut
- Powered by a keydrive
- Motorsold separately

Item	Product Description	UOM
53089	Power Handle™ Rotary Float, 17-1/2′Length,90Degree, 1"DiamondDisc	Each
53090	Power Handle™ Rotary Float, 7–1/2"Length,Cylindrical Bur	Each
53091	Power Handle™ Rotary Float, 10"Length, Cylindrical Bur	Each
53092	Power Handle™ Rotary Float, 17-1/2'Length,Cylindrical Bur	Each
53093	Cylindrical Bur, 5/81Diameter, 2"Length, Carbide	Each
53094	Cylindrical Bur, 5/8'Diameter, 2"Length, Diamond	Each



Tri-TipGrinder Kits Poental Products

The Tri-Tip Grinder Kit is a unique grinder kit used to smoothout the molars, incisors and canine teeth in all areas of the dental arcade from miniature to draft horses. The knurled ultra slimprofile quard measures 18" long. The kit indudes your choice of interchangeable protective shrouds measuring 3/8', 1/2'or 5/8' diameter and a removable knurled handle measuring 6-1/4' long. Also offered is a choice of a 6' long cylindrical diamondchip bur 3/8', 1/2" or 5/8" in diameter and a 1" handpiece. The grinder is powered by a key drive external motorand may be used with any power unit. Dreme® or Foredom[™]. Additional shrouds and burs are available separately. Motors also sold separately.

United States: 1.800.447.0687

ltem	Product Description	UOM
	Tri-TipCrinder Kit, 1/2'Shroud, 18'Length	Each
	Tri-TipCrinder Kit, 3/8'Shroud, 18'Length	Each
	Tri-TipCrinder Kit, 5/8'Shroud, 18'Length	Each
53387	Cylindrical Bur, 1/2'Diameter,6"Length, Split Blue Diamond	Each
	Cylindrical Bur, 3/8'Diameter,6"Length, Split Blue Diamond	Each
53389	Cylindrical Bur, 5/8'Diameter,6"Length, Split Blue Diamond	Each




DEWALT is a registered trademark of DeWALT Industrial Tool Company. DREMEL is a registered trademark of Robert Bosch Tool Corporation.







OuickHoat™ Kit JERRY

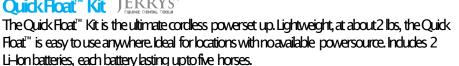

Kit indudes:

- Quick Float[™] Motor
- QuidkFloat™ Wandwith choice of DiamondBur
- 2 Recharcable Li-Ion Batteries
- Wall ChargerUnit
- Lubricant
- LubeGun
- **DiskPliers**
- T-Wrench
- Convenient Carrying Case

ltem	Product Description	UOM
58100	Quidk Float™ Kit	Each



EQ Power Float[™]



The EQ Power Float™ measures 19 inches in length and features a low profile 90 degree head. The EQ Power Float™ is significantly lighter than most on the markettoday, weighing under 1.5 lbs. The EQ Power Floats are built in-housewith high quality aluminum and stainless steel. Available with choice of burand choice of Key Drive or Square Drive.

Item	Product Description	UOM
58102	EQ Power Float™ -Key Drive	Each
58103	EQ Power Float™ -Square Drive	Each



EQMiniPower Hoat JERRYS



The EQ Mini Power Float™ measures 8.5 inches in length and features a low profile 90 degree head. Small and easy to maneuver, the EQ Mini Power Float™ is ideal for work on incisors and canines. The EQ Power Floats are built in-housewith highquality aluminumandstainless steel. Available with choice of burand choice of Key Drive or Square Drive.

ltem	Product Description	UOM
58104	EQ Mini Power Float™ -Key Drive	Each
58106	EQ Mini Power Float™ -Square Drive	Each



MightyEZ Hoat™ Kits JERRY



The MightyEZ Float™ is a lightweight, battery-poweredunit. Featuring an On/Offswitch, variable speed control and a built in dutch, the MightyEZ Float™ offers all of the standard features in a convenient, cordless unit.

Kit indudes:

- 2 Rechargeable Li-lon Batteries
- Wall ChargerUnit
- 43" or 66" Key Drive or Square Drive Complete Cable Assembly
- Convenient Carrying Case

ltem	Product Description	UOM
58110	MightyEZ Float™ Kit -Key Drive, 43"	Each
58111	MightyEZ Float™ Kit -Square Drive, 43"	Each
58108	MightyEZ Float™ Kit -Key Drive, 66"	Each
58109	MightyEZ Float™ Kit -Square Drive, 66'	Each









The Makita® Motor is a powerful, reliable unit requiring minimal maintenance.

- 43" or 66" Key Drive or Square Drive Complete Cable Assembly
- Foredom™FootPedal, Plastic

Item	Product Description	UOM
58049	Makita® Motor	Each
580498	Makita® Motor-Short	Each





ForecomTM Motors Poental Products

- Key drive powershaft is compatible with the discontinued Dremel® products
- Variable speed controlled by foot pedal
- Heavy duty, high power unit
- Delivers moretorque with complete control throughout the entire speed range allowing for greater control of the bur

ltem	Product Description	UOM
58058	Foredom™TX-SXR 1/3hpMotorwithFoot Pedal -Key Drive	Each
58059	Foredom™SR-SCT 1/6hpMotorwithFoot Pedal -Square Drive	Each
58075	Foredom™63"Inner Cable, Square Drive	Each
58076	Foredom™63"Outer Cable, Square Drive	Each
58077	Foredom™39"Inner Cable, Key Drive	Each
58079	Foredom™39'Outer Cable, Key Drive	Each
58078	Foredom™66"Inner Cable, Key Drive -Non-conductive	Each
58080	Foredom™66"Outer Cable, Key Drive	Each
58073	Foredom™45-3/4³hnerCable, Square Drive -Non-conductive	Each
58074	Foredom™45-3/4"OuterCable, Square Drive -Non-conductive	Each



Power Motor & Accessories Power Motor & Accessories Power Motor & Accessories

- Lowcost, high value power motorcomparable in durability to Foredom™
- Powerful, compactand lightweight
- Available in either 110-Voltor220-Volt


Item	Product Description	UOM
53096	Power Motorwith Foot Pedal	Each

EZY-Drive™ Motor System ✓ Preference Ary Dental Products.

A compact, lightweight all-in-one battery-powered float system for equine dentistry. Features:

- Quiet motorunit that provides adequate torque
- Independent systempower On/Offswitch
- Independent variable speed control dial
- Operating range of 1,200-9,500 RPM
- Li-Ion battery, 25.9-Volt, 3.6 amp-hours
- Battery fully charges in 1.5 hours
- Battery runtimelasts an average of 10 horses
- Keyed drive -use with any key-driven handpiece

Item	Product Description	UOM
53155	EEZY-Drive™ MotorSystem-Base Model	Each
53156	EEZY-Drive™ MotorKit	Each
53157	EEZY-Drive™ Harness	Each

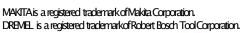









53155



Equine Dental Systems Polyfloat & Accessories Polyfloat & Polyfloat & Dental Products...

. , , , , , , , , , , , , , , , , , , ,		
Item	Product Description	UOM
53163	EDS Polyfloat Straight, Key Drive, Short -15 cmShaft Length (without bur)	Each
53164	EDS Polyfloat Straight, Key Drive, Standard –25 cmShaft Length (without bur)	Each
53165	EDS Polyfloat: Straight, Key Drive, Long-37cmShaft Length (without bur)	Each
53166	EDS Polyfloat Curved, Key Drive, Standard -25 cmShaft Length (without bur)	Each
53167	EDS Polyfloat Curved, Key Drive, Long-37 cm Shaft Length (without bur)	Each




Polyfloat Burs Polyfloat Burs Polyfloat Burs

ltem	Product Description	UOM
53168	Polyfloat: Applecore Bur, Canine –18 mm Diameter	Each
53169	Polyfloat: Applecore Bur -14 mmDiameter	Each
53170	Polyfloat Applecore Bur -18 mmDiameter	Each
53171	Polyfloat: Applecore Bur -22 mmDiameter	Each
53172	Polyfloat: Ocdusal Bur -17 mmDiameter	Each
53173	Polyfloat Cylindrical Bur -14mmLength	Each
53174	Polyfloat Cylindrical Bur -27mmLength	Each
53175	Polyfloat Versatile Cylindrical Bur –20 mmLength	Each
53176	Polyfloat Versatile Cylindrical Bur -25 mmLength	Each
53177	Polyfloat Conical Bur -10 mmLength	Each
53178	Polyfloat Conical Bur -15 mmLength	Each



Equine Dental Systems Flexicisc Poental Products...

	ι
e, Short -15 cmShaft Length	E
e, Standard –25 cmShaft Length	E
e Long-37cmShaft Length	E

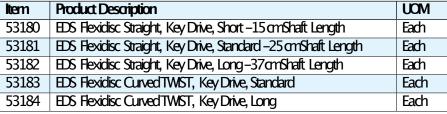



Handpiece Drils Poental Products

Description	
Dill	
ılary Drill	
1	





ltem	Product Description	UOM
53192	Extension, Long	Each
53193	Holster, Right	Each
53194	Holster, Left	Each







Equine Periodontal System& Accessories Poental Products...

The Equine Periodontal Systemis used for the diagnosis and treatment of periodontal disease and endodontic procedures in horses. Systemis complete and ready to use -the user only needs to supplyatankofnitrogen.

High Speed Handpiece with an equine extension

- Air Water Handpiece with both short and long tips
- Air Abrasion Handpiece
- Vacuumhandpiece
- Can be operated using Nitrogen or compressed air
- Nitrogen regulator may be purchased separately
- Mobile cart may be purchased separately

Item	Product Description	UOM
53062	Equine Periodontal System	Each
53060	Equine Periodontal Systemwith Cart	Each
53061	Equine Periodontal Cart	Each

Air Abrasion Handpiece P Dental Products

Item	Product Description	UOM
53068	Air Abrasion Handpiece	Each

High Speed, High Torque Handpiece Poental Products...

Item	Product Description	UOM
53130	High Speed, High Torque Handpiece -350,000 RPM	Each

High Speed Extension Products

Item	Product Description	UOM
53069	High Speed Extension -Midwest 4 Hole Coupling	Each



•		
ltem	Product Description	UOM
53131	Replacement Turbine for High Speed, High Torque Handpiece	Each

Prophy Jet Powder Prophy Jet Powder Powder Prophy Jet Powder

Item	Product Description	UOM
53749	Prophy Jet Powder, 600 grams	Each

Tapered Fissure #702L, Surgical, FG, 25mm-Carbide Dental Products ac

ltem	Product Description	UOM
54812	Tapered Fissure #702L, Surgical, FG, 25 mm-Carbide	5-Pack

Round Tip#8, FG, Surgical, 25mm-Carbicle Poental Products...

Item	Product Description	UOM
53042	Round Tip#8, FG, Surgical, 25mm-Carbide	5-Pack

3-WaySyringe with Short Tip Pental Products.


Item	Product Description	MOU
53065	3-WaySyringe with Short Tip	Each









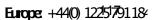



Extended Syringe Tip VETERINARY Dental Products...

ltem	Product Description	UOM
53067	Extended Syringe Tip -Tip Only	Each

Bio Lube[®] Lubricant Pen Poental Products...

Item	Product Description	UOM
54882	Bio Lube®Lubricant Pen, 25oz	Each



Gunther Speculum

Item	Product Description	LOM
53612	Gunther Speculum	Each





Pacific Edge™ MouthSpeculum ♥ Dental Products...

Item	Product Description	UOM
53610	Pacific Edge™ MouthSpeculumwith Biothane Straps, Polished Finish	Each



Harton's MouthWedge Pental Products

 Unique design allows several teeth to rest on the wedge at a time, reducing the danger of gradking

ltem	Product Description	UOM
53689	Harlton's MouthWedge	Each



Equine Cheek Retractor Pental Products

This Equine Cheek Retractor, or spreader, is designed to gently pull the cheeks away from the cheek teeth so occlusion can be evaluated on a dosed mouth. It is also useful for providing more room when using power dentistry equipment. The stainless steel design features two offsetting size retractors, one large and one small, that can be easily alternated by simply turning the retractor around.

Item	Product Description	UOM
53615	Equine Cheek Retractor	Each



##EZY-View Dental Carnera & Accessories • VETERINARY Dental Products...

- Crisp, dear image
- LED lights (6 levels)
- 360-degreerotation of camera shaft
- Portable unit
- Detachable, wireless LCD monitor
- Micro SD Card included for image /video storage

Item	Product Description	UOM	
52500	EEZY-View™ Base Unit &LCD Monitorwith Carrying Case	Each	
52505	EZY-View™ Camera, 10mmODx430mmLength, Rigid, 90-DegreeView	Each	
52510	EZY-View™ Camera, 10mmODx430mmLength, Semi-Rigid, 90-DegreeView	Each	
52520	EEZY-View™ Camera, 6mmODx 1000mmLength, Semi-Rigid, Front View	Each	
52515	EEZY-View™ LCD Monitor	Each	
52517	EEZY-View™ Monitor Replacement Charging Cord	Each	
52519	EEZY-View™ Camera, 32GB MicroSD Card with Adapter	Each	
ThoEdo	The February Water Pick Kits & Accessories The Edge		





The Edge Equine™ Water Pick Kits & Accessories

- Ideal for dislodging impacted food from diasternata.
- Madefromacetyl, stainless steel and brass

Item	Product Description	UOM
53762	The Edge Equine™ Quick ChangeWater Pick & Flushing Wand Kit with Bottle	Each
53763	The Edge Equine™ Quick Change Water Pick & Flushing Wand with Tubing Assembly	Each



Wide Spray Pattern





53762

Equine Radiography Top VETERINARY Dental Products

ltem	Product Description	UOM
53427	Equine RadiographyTop(fits Standard & Cradle Style Base)	Each



Equine HeadRest, Standard Style Poental Products.

- Constructed of lightweight aluminum
- Features a tip resistant base, twist-proofsquare shankand height adjustable lock
- Powdercoated forease in dean up
- Quidkly disassembled for transportation or storage
- Traditional roundstyle resting pad
- Compactandeasy to store

Ite m	Product Description	UOM
53430	Equine Head Rest & Base, Standard Style	Each
53437	Replacement Cover & Pad, Standard Style	Each



Equine HeadRest, Cracle Style Products

- Constructed of lightweight aluminum
- Features a tip resistant base, twist-proofsquare shankand height adjustable lock
- Powder coated for ease in dean up
- Quidkly disassembled for transportation or storage
- Large paddedplatformprovides a large resting area for horse's head
- Platform with gradually sloping sides holds the sedated horse more securely and places less strain on the mandible
- Compactandeasy to store

ltem	Product Description	UOM
53433	Equine Head Rest & Base, Cradle Style	Each
53434	Replacement Cover & Pad, Cradle Style	Each



QuikLift Equine HeadRest, Standard Style Poental Products.

Self-elevating head rest makes holding a sedated horse head, while adjusting the head rest height and angle, a one-handedoperation.

- Constructed of lightweight aluminum
- Features a tip resistant base, twist-proofsquare shankand height adjustable look
- Powdercoated forease in dean up
- Quidkly disassembled for transportation or storage
- Traditional roundstyle resting pad
- Compactandeasy to store
- Self rising

Item	Product Description	UOM
53430QL	Quilk Lift Equine Head Rest, Standard Style	Each
53437	Replacement Cover & Pad, Standard Style	Each











QuikLift Equine HeadRest, Cradle Style Poental Products

Self-elevating head rest makes holding a sedated horse head, while adjusting the head rest height and angle, a one-handedoperation.

- Constructed of lightweight aluminum
- Features a tip resistant base, twist-proofsquare shankand height adjustable look
- Powder coated for ease in dean up
- Quidkly disassembled for transportation or storage
- Traditional roundstyle resting pad
- Compactandeasy to store
- Self rising

	ltem	Product Description	UOM
	53433QL	Quik Lift Equine Head Rest, Cradle Style	Each
Ì	53434	Replacement Cover & Pad, Cradle Style	Each



QuikLift Conversion Post Poental Products

- Constructed of lightweight aluminum
- Powder coated for ease in dean up
- Self rising
- Converts an original head rest to a Quik Lift Head Rest

ltem	Product Description	UOM
53438	QuikLift ConversionPost (only fits Standard & Cradle Style Base)	Each



Equine Cactaver Head Holder Pental Products

Sturdy holderfor cadaver heads used for educational purposes.

Item	Product Description	UOM
53401	Equine Cadaver Head Holder(fits Standard & Cradle Style Base)	Each
53402	Equine Cadaver Head Holder & Base	Each

Equine HeadRest, SuspendedStyle Poental Products... PATENTED

- Provides excellent and safe support for the sedated horse's head during all dental procedures, opthalmic exams or danial surgeries
- Provides adjustable cranial elevation that is self-locking, secure and stable, yet allows formild swaying and is comfortable and safe for the patient
- Designed to be suspended from a beamor rafter using a chain or rope or suspended from stods.
- Camdeat rope adjuster allows quidk variable height adjustments using a self-locking quidk release mechanism
- Indudes a high quality foampad, vinyl cover and hardware.

Item	Product Description	UOM
53431	Equine Head Rest, Suspended Style	Each
53436	Replacement Cover & Pad, Suspended Style	Each



Dental Halter System JERRY'S

The Dental Halter Systeminducies a dental halter and fiddle block along with a pulley and rope.

	,	<u> </u>
ltem	Product Description	UOM
51802	Dental Halter System	Each





LaGrangeS-ShapedScissors -Senated



Item	Product Description	UOM
54580	LaGrange S-Shaped Scissors -Serrated	Each



DeanScissors



ltem	Product Description	UOM
54581	Dean Scissors, 7'Length	Each



Iris Scissors -SuperCut:CurvedSerrated



Item	Product Description	UOM
54582	Iris Saissors -SuperCut Curved Serrated, 4.75"Length	Each



Iris Scissors -SuperCutStraight Serrated



Item	Product Description	UOM
54584	Iris Scissors -SuperCut Straight Serrated, 4.75"Length	Each



Olsen Hegar Needle Holder with Scissors



	"Meal instead"	
Item	Product Description	UOM
54583	Olsen Hegar Needle Holderwith Sassors, 5.5" Length	Each



MayoScissors, Straight 6.75"



Mayo Scissors are a versatile instrumentused in multiple procedures and settings. These scissors are primarily used in procedures to cut heavier tissue. They are available with either straight or curved blades.

Г	- -	Durch at Dana dadan	1100.4
ı	ltem	Product Description	UOM
ſ	75410SG	Mayo Scissors, Straight 6.75"Length	Each



Olsen Hegar Needle Holder, Tungsten Carbide USI



Olsen Hegar Needle Holder is an "all-in-one" needle holding instrument that features cross-serated tips to hold the suture needle in addition to a built-inscissors to then cut the suture. This needle holder features tungsten carbide tips, gold finger ring handles and are madewith high quality Germanstainless steel.

ltem	Product Description	UOM
75616G	Olsen Hegar Needle Holder; Tungsten Carbide	Each



SHORTY[™] Sutures

5.10.11.										
PV		POLYVINYLIDENE FLUORIDE, NON-AE								
ABSORPTION	ON PROF	lle .	Exceptional Knot Security, knots fla	atten and lo	dk against themselves, unwetted by blood, inert Non	-Absorbal	ole ole			
ORDER CODE	ODER CODE USP GALICE EP NEEDLE TYPE NEEDLE SIZE		NEEDLE	SUTURE LENGTH	OTY/ BOX					
PV20425D	2–0	3	3/8 Reverse Cutting $ abla$	25mm		45cm	36			
PV30425D	3–0	2	3/8 Reverse Cutting ▽	25mm		45cm	36			

Double forward throw

(Surgeons Throw)



Single reverse throw

Creating a square knot, with which to adjust tension.



Single forward throw

(Locking the knot)



Final result

9 interrupted sutures, with long tails, from a 45cm length.



Attention A Veterinarians...



SHORTY™

Sutures rebased

Ever just needed a short length high vis suture for a small incision?

Look no further than SHORTY™.

Low Cost – No Waste.

Special introductory offer, for a limited time!

Innovation/VetsTalk,SlverClicle Listens

Fluoro Pink Polypropylene Sutures

RR		ENT PINK MONOFILAMENT, NON-/	ABSORB	ABLE			
ABSORPTION PROFILE		LE .	Unwetted by blood, inert, knots fla		k against themselves for improved knot security NO	N-ABSORB	ABLE
ORDER CODE	USP GAUGE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	OTY/ BOX
RR2075D	2-0	3	1/2 Cirde Taper Point ⊙	26.5mm	\	70cm	36
RR0075D	0	3.5	1/2 Gide Tapa Foilite	20.311111		70cm	36
RR50422N	5–0	1	3/8 Circle Reverse Cutting $ abla$	16mm		45cm	12
RR4046N	4-0	1.5	2/9 Girdo Dorrano Gittings	10,000		45cm	12
RR4046D	4-0	1.5	3/8 Cirde Reverse Cutting ▽	18mm		45cm	36
RR3077N	3-0	2	3/8 Girde Reverse Cutting ▽			70am	12
RR3077D	3-0	2		24mm		70am	36
RRCLK3047D [†]	3–0	2				45cm	36
RR20710D	2–0	3	3/8 Cirde Reverse Cutting ▽			70am	36
RRCLK20410D [†]	2-0	3		30mm		45cm	36
RRCLK00310D	0	3.5				45cm	36
RR00113D	0	3.5	1/2 Circle Reverse Cutting $ abla$	37mm		100cm	36
RR00130D	0	3.5				100cm	36
RR01130D	1	4				100cm	36
RR02130D	2	5	2/9 Guda Darama Guttinav√	F0		100cm	36
RRCLK00430D [†]	0	3.5	3/8 Cirde Reverse Cutting ▽	50mm		45cm	36
RRCLK01430D [†]	1	4				45cm	36
RRCLK02430D [†]	2	5				45cm	36
RR30760D	3-0	2	Stanislat C. Hings			70am	36
RR20760D	2-0	3	Straight Cutting ▽	60mm	_	70am	36

† Indicated for Caslick Surgery and small incisions

Examples of high visibility sutures making for easier suture placement and removal.








Polypropylene Sutures

PP		BLUE MONOFILAMENT, NON-ABSORBABLE								
ABSORPTION F	ROFILE		Unwetted by blood, inert, knots flat	ten and loc	k against themselves for improved knot security No r	n-Absorbal	ble			
ORDER CODE	USP GAUCE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	QTY/ BOX			
PP3077D	3–0	2	3/8 Circle Reverse Cutting $ abla$	24mm		70am	36			
PP20710D	2-0	3	3/8 Circle Reverse Cutting ∇	30mm		70cm	36			

Braided Polyester Sutures

TieBack Equine

BP		BRAIDED POLYESTER, GREEN SILICONE COATED,NON-ABSORBABLE								
ABSORPTION PROFILE			High strength braid Non-Absorba	able						
ORDER CODE	USP GAUCE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	OTY/ BOX			
BP05737N	5	7	1/2 Circle Reverse Cutting $ abla$	37mm		75am	12			
BP05738N	5	7	1/2 Cirde Taper Point ⊙	37mm		75am	12			

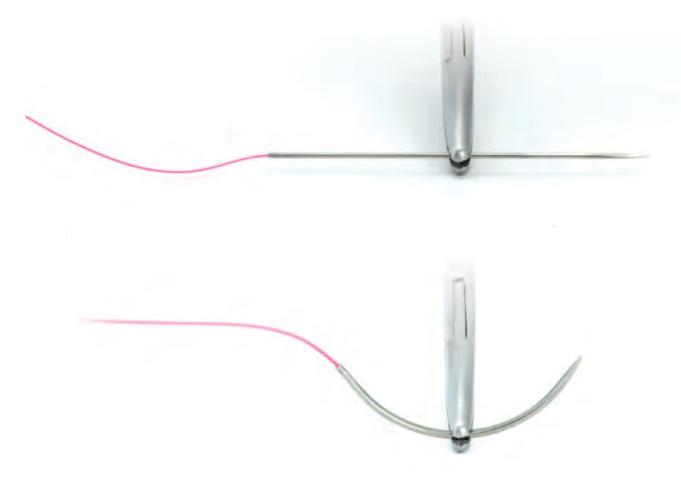

Nylon Sutures

N					BLACK MONOFILAMENT, NON-ABSO	RB	ABLE
ABSORPTION	PROFILE		High strength retention Non-Abs	orbable			
ORDER CODE	USP GAUCE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE SUTT	IRE ITH	OTY/ BOX
N60416N	6-0	0.7	3/8 Cirde Reverse Cutting $ abla$	11mm	450	m	12
N50417N	5–0	1	3/8 Cirde Reverse Cutting $ abla$	12mm	450	m	12
N50422N	5-0	1	2/9 Circle Day area Cuttings	16mm	450	m	12
N50422D	5-0	1	3/8 Circle Reverse Cutting ▽	TOULIU	450	m	36
N4046N	4-0	1.5	3/05/ds Paraga 6 #ins	10	450	m	12
N4046D	4-0	1.5	3/8 Circle Reverse Cutting ▽	18mm	450	m	36
N4047N	4-0	1.5			450	m	12
N4077D	4-0	1.5	3/8Cirde Reverse Cutting ▽		700	m	36
N3047N	3-0	2		24	450	m	12
N3047D	3-0	2		24mm	450	m	36
N3077D	3-0	2			700	m	36
N3077N	3-0	2			700	m	12
N20710D	2-0	3			700	m	36
N00110D	0	3.5	3/8Cirde Reverse Cutting $ abla$	30mm	100	m	36
N01140D	1	4			100	m	36
N01130D	1	4	3/8 Girde Reverse Cutting $ abla$	50mm	100	m	36
N30760D	3-0	2	S : 1. S :: \(\tau \)	60	700	m	36
N20760D	2-0	3	Straight Cutting ▽	60mm	700	m	36
N02125D	2	5	3/8 Girde Reverse Cutting ▽	<i>77</i> mm	100	m	36

Fluoromid[™] Sutures

SMF			FLUORESCENT YE	LOW N	MULTIFILAMENT POLYAMIDE, NON-	ABSORE	ABLE			
ABSORPTIO	ON PROF	le .	Easy handling, ultra high visibility Non-Absorbable							
ORDER CODE	USP CAUCE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	OTY/ BOX			
SMF3077D	3-0	2	3/8 Cirde Reverse Cutting $ abla$	24mm		70am	36			
SMF20710D	2-0	3	- 3/8 Circle Reverse Cutting $ abla$	30mm		70am	36			
SMF00110D	0	3.5		SOTTI		100cm	36			
SMF00113D	0	3.5	1/2 Circle Reverse Cutting $ abla$	37mm		100cm	36			
SMF01130D	1	4	3/8 Circle Reverse Cutting $ abla$	50mm		100cm	36			
SMF20760D	2-0	3	6 : 1.6 :: 5	60mm		70am	36			
SMF30760D	3-0	2	Straight Cutting $ abla$	WIIII		70am	36			

White Fluoromid[™] Cassettes

White Fluromid [™]									
ABSORPTIONPROFILE Easy handling Non-Absorbable									
CODE	USP GAUGE	P	LENGTH	QTY					
SMCAS40100	4-0	1.5	100	1					
SMCAS30100	3–0	2	100	1					
SMCAS20100	2-0	3	100	1					
SMCAS00100	0	3.5	100	1					
SMCAS01050	1	4	50	1					
SMCAS02050	2	5	50	1					

127
PAGE 7 | SILVERGLIDE CATALOGUE | JULY 2019

Polydioxanone Sutures

	DI IDDI E MONOEII AMENITAPSODPADI E								
PDO					PURPLE MONOFILAMENT, ABSOI	RB/	ABLE		
ABSORPTIO		LE	50-80% retained tensile strength		Essentially complete between 6-9 months.				
ORDER CODE	USP GAUGE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE SUTUL LENG	H	QTY/ BOX		
PO50724D	5-0	1	1/2 Cirde Taper Point ⊙	16mm	70cr	n	36		
PO4074D	4-0	1.5	1/2 Cirde Taper Point ⊙	22mm	700	n	36		
PO3075D	3-0	2			70a	n	36		
PO2075D	2-0	3	1/2 Cirde Taper Point ①	26.5mm	70cr	n	36		
PO0015D	0	3.5			1000	n	36		
PO2078D	2-0	3			70cr	n	36		
PO0118D	1	4	1/2 Cirde Taper Point ⊙	36mm	1000	n	36		
PO0018D	0	3.5			1000	n	36		
PO011339N	1	4			1300	n	12		
PO021339N	2	5	1/2 Cirde Taper Point ⊙	65mm	1300	n	12		
PO03139N	3	6			1000	n	12		
PO4076D	4-0	1.5	3/8 Cirde Reverse Cutting ▽	18mm	7000	n	36		
PO3077D	3-0	2	3/8 Cirde Reverse Cutting $ abla$	24mm	7000	n	36		
PO20726D	2-0	3	3/8 Girde Reverse Cutting $ abla$	26mm	7000	n	36		
PO20713D	2-0	3	1/2 Girde Reverse Cutting $ abla$	37mm	70cr	n	36		
POCLK20410D [†]	2-0	3	3/8 Girde Reverse Cutting $ abla$	40mm	450	n	36		
PO00114D	0	3.5			1000	n	36		
PO01114D	1	4	1/2 Circle Reverse Cutting $ abla$	40mm	1000	n	36		
PO02114D	2	5			1000	n	36		
PO01130D	1	4			1000	n	36		
PO02130D	2	5			1000	n	36		
POCLK00430D [†]	0	3.5	- 3/8Cirde Reverse Cutting▽	50mm	450	n	36		
POCLK01430D [†]	1	4	3/0 Glac Neverse Culling V	JUITIII	450	n	36		
POCLK02430D [†]	2	5			450	1	36		
PO001130D	0	3.5			1000	n	36		

† Indicated for Caslick Surgery and small incisions

Polyglycolic Acid Sutures

PGA					BRAIDED VIOLET, AB	SORB/	ABLE.
ABSORPTIO	ON PROF	LE .	65% retained tensile strength at	14 davs. Es	sentially complete between 3–6 months.		
ORDER CODE	USP	EP	NEEDLE TYPE	NEEDLE SIZE		SUTURE LENGTH	QTY/ BOX
PG4074D	4-0	1.5	1/2 Cirde Taper Point ⊙	22mm		70cm	36
PG3074D	3–0	2	1/2 circle raper rollic	2211111		70am	36
PG3075D	3-0	2			_	70cm	36
PG2075D	2–0	3	1/2 Cirde Taper Point ⊙	26.5mm		70cm	36
PG0075D	0	3.5				70cm	36
PG0115D	1	4				100cm	36
PG0018D	0	3.5	1/2 Cirde Taper Point ⊙	36mm		100cm	36
PG0118D	1	4				100am	36
PG00140D	0	3.5	1/2 Cirde Taper Point ⊙	40mm		100cm	36
PG01140D	1	4	-/- C			100cm	36
PG011339N	1	4				130cm	12
PG021339N	2	5	1/2 Cirde Taper Point ⊙	65mm		130cm	12
PG031339N	3	6	1/2 Circle Taper Point.	OSMIN		130cm	12
PG03139N	3	6				100cm	12
PG031390N	3	6	3/8 Cirde Taper Point ⊙	90mm		130cm	12
PG60417N	6-0	0.7	2/05-d- Parent C #ins \	12		45cm	12
PG50417N	5-0	1	3/8Cirde Reverse Cutting ▽	12mm		45cm	12
PG4076D	4-0	1.5	3/8 Cirde Reverse Cutting ▽	18mm		70am	36
PG3077D	3–0	2	2/9 Giralo Der mine Carrier	24-2		70cm	36
PGCLK3047D [†]	3-0	2	3/8 Circle Reverse Cutting ▽	24mm		45cm	36
PG20721D	2-0	3	1/2 Girde Reverse Cutting $ abla$	26mm		70am	36
PG00110D	0	3.5				100cm	36
PGCLK20410D [†]	2-0	3	3/8 Circle Reverse Cutting $ abla$	30mm		45cm	36
PGCLK00410D [†]	0	3.5				45cm	36
PGCLK01430D [†]	1	4	2/9 Circle Deverse Cotting	FO		45cm	36
PGCLK02430D [†]	2	5	3/8 Circle Reverse Cutting $ abla$	50mm		45cm	36

† Indicated for Caslick Surgery and small incisions

Polyglycolic Acid Fast AbsorbingSutures

Dental

PGF		///		////	BRAIDED, UNDYED,	ABSORBA	ABLE
ABSORPTIO	entially complete at 14 days.						
ORDER CODE	USP GAUCE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	QTY/ BOX
PGF50413N	5–0	1	1/2 Cirde Reverse Cutting $ abla$	13mm	\bigcup	45cm	12
PGF40422N	4-0	1.5	3/8 Cirde Reverse Cutting $ abla$	16mm		45cm	12
PGF3046N	3-0	2	3/8 Cirde Reverse Cutting $ abla$	18mm		45cm	12



Polyglactin Sutures

PGLA					BRAIDED WOLET, A	BSORB	ABLE
ABSORPTION PROFILE 74% retained tensile strength at 14 days. Essentially complete at 63 days.							
ORDER CODE	USP GAUGE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	QTY/ BOX
PGLA03939N	3	6	1/2 Cirde Taper Point ⊙	65mm		90cm	12



Glycolide & eCaprolactone Sutures

CCL	PURPLE MONOFILAMENT, ABSORBABLE							
ABSORPTION	ON PROF	ILE .	60% retained tensile strength at 2	7 days. Esse	entially complete between 2-4months.			
ORDER CODE	USP GAUGE	Ð	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	QTY/ BOX	
EC50450N	5–0	1	1/2 Girde Taper Point ⊙	11mm	\bigcup	45cm	12	
EC40416D	4-0	1.5	1/2 Cirde Taper Point ⊙	16mm		45cm	36	
EC50428N	5-0	1	1/2 Cirde Taper Point ⊙	18mm	\	45cm	12	
EC40728D	4-0	1.5	1/2 Glue Tapa Politic	1011111		70cm	36	
EC4074D	4-0	1.5				70cm	36	
EC3074D	3-0	2	1/2 Cirde Taper Point ⊙	22mm		70cm	36	
EC2074D	2-0	3				70cm	36	
EC3075D	3-0	2	1/2 Girda Tanar Paint	26.5mm	\ 1	70am	36	
EC2075D	2-0	3	1/2 Cirde Taper Point ⊙	20.311111		70am	36	
EC0018D	0	3.5	1/2 Girda Tanar Paint	26,000	\ /	100am	36	
EC0118D	1	4	1/2 Cirde Taper Point ⊙	36mm		100am	36	
EC011339N*	1	4	1/2 Cirde Taper Point ⊙	65mm		130cm	12	
EC50416N	5–0	1	3/8 Cirde Reverse Cutting $ abla$	16mm		45am	12	
EC4076D	4-0	1.5	3/8 Cirde Reverse Cutting $ abla$	18mm		70am	36	
EC3077D	3-0	2	3/8 Cirde Reverse Cutting ▽	24mm		70cm	36	
EC20726D	2-0	3	3/8 Cirde Reverse Cutting ▽	26mm		70cm	36	
EC20713D	2-0	3	1/2 Circle Por over Cuttings √	27	\	70cm	36	
EC00113D	0	3.5	1/2 Circle Reverse Cutting $ abla$	37mm		100cm	36	

^{*}Bovine uterus dosure only

Sutures Non-Needled Absorbable

	Polyglyd	colic Aci	d					
ABSORPTION PROFILE		65 % retained tensile strength at 14 days . Essentially complete between 3–6months						
CODE	USP GAUGE	P	LENGTH	QTY/ BOX				
PGNN3045D	3-0	2	45	36				
PGNN2045D	2-0	3	45	36				
PGNN0045D	0	3.5	45	36				
PGNN0145D	1	4	45	36				
PGNN0245D	2	5	45	36				
PGNN30100D	3-0	2	100	36				
PGNN20100D	2-0	3	100	36				
PGNN00100D	0	3.5	100	36				
PGNN01100D	1	4	100	36				
PGNN02100D	2	5	100	36				

	Glycolid	e & eCa	р						
ABSORPTION PROFILE		60% retained tensile strength at 7 days Essentially complete between 2-4 months							
CODE	USP GAUGE	P	LENGTH	36 36 36 36 36 36 36					
ECNN3045D	3-0	2	45	36					
ECNN2045D	2-0	3	45	36					
ECNN0045D	0	3.5	45	36					
ECNN0145D	1	4	45	36					
ECNN0245D	2	5	45	36					
ECNN30100D	3-0	2	100	36					
ECNN20100D	2-0	3	100	36					
ECNN00100D	0	3.5	100	36					
ECNN01100D	1	4	100	36					
ECNN02100D	2	5	100	36					

	Polydioxanone								
ABSORPTION PROFILE	42 days. Es	50-80% retained tensile strength at 42 days. Essentially complete between 6-9 months							
CODE	USP GAUGE	P	LENGTH	OTY/ BOX					
PONN3045D	3-0	2	45	36					
PONN2045D	2-0	3	45	36					
PONN0045D	0	3.5	45	36					
PONN0145D	1	4	45	36					
PONN0245D	2	5	45	36					
PONN30100D	3-0	2	100	36					
PONN20100D	2-0	3	100	36					
PONN00100D	0	3.5	100	36					
PONN01100D	1	4	100	36					
PONN02100D	2	5	100	36					

Reusable Staplers

CODE	DESCRIPTION	QUANTITY
SRH-1	Pre-Loaded SRH Handle	1
SRC35W	Cartridges	1
SRE	Staple Extractor	1



Bovine Surgery

Fluoromid[™] Sutures

SMF		FLUORESCENT YELLOW OR WHITE, NON-ABSORBABLE								
ABSORPTION PROFILE Easy handling, ultra high visibi			Easy handling, ultra high visibility	Non-Abso	rbable					
ORDER CODE	USP GAUGE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	QTY/ BOX			
SMF031360N	3	6	Straight Cutting $ abla$	65mm	_	130am	12			
SMF031375N	3	6	$1/2$ Circle Conventional Cutting / Straight Shank (Ski) \triangle	75mm		130am	12			

Glycolide & eCaprolactone Sutures

Uterus Closure

GCL		PURPLE MONOFILAMENT, ABSORBABLE								
ABSORPTIO	ON PROF	ile.	entially complete between 2-4months.							
ORDER CODE	USP GAUGE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	OTY/ BOX			
EC011339N	1	4	1/2 Circle Taper Point ⊙	65mm		130am	12			

Polyglycolic Acid Sutures


PGA					BRAIDED VIOLET, A	BSORB	ABLE
ABSORPTIO	sentially complete between 3-6months.						
ORDER CODE	USP GAUCE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	OTY/ BOX
PG011339N	1	4	1/2 Cirde Taper Point ⊙	65mm		130cm	12
PG021339N	2	5	1/2 Cirde Taper Point ⊙	65mm		130cm	12
PG031339N	3	6	1/2 Cirde Taper Point ⊙	65mm		130cm	12

Polydioxanone Sutures


PDO					PURPLE MONOFILAMENT,A	BSORB	ABLE
ABSORPTION PROFILE 50-80% retained tensile strength at 42 days. Essentially complete between 6-9 months					Essentially complete between 6-9months.		
ORDER CODE	USP GAUGE	EP	NEEDLE TYPE	NEEDLE SIZE	NEEDLE	SUTURE LENGTH	QTY/ BOX
PO021339N	2	5	1/2 Cirde Taper Point ⊙	65mm		130cm	12

Needle Chart

Needle Shapes

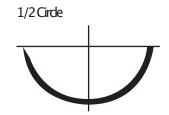

Straight Cutting

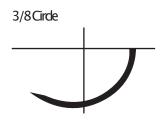

1/2 Circle Conventional Cutting / Straight Shank (Ski)

Needle Types



Round bodied taper point needle \odot



Conventional cutting needle \triangle



Reverse cutting needle ∇

Circle Reverse Cutting

1/2Circle Reverse Cutting ∇

Needle	Size	SilverGlide™	Ethicon ™ (similar)	Syneture [™]	B.Braun [™] (similar)
	\bigtriangledown 13mm	C13	-	-	-
	▽ 26mm	C21	CP-2	CS-10	HS26
	▽ 37mm	C37	CP-1	GS-11	HS37
	▽ 40mm	C14	œ	CS-12	HS40

3/8Circle Reverse Cutting ∇

Needle	Size	SilverGlide™	Ethicon [™] (similar)	Syneture™	B.Braun [™] (similar)
	\bigtriangledown 11mm	C16	C-2	C-1	DS12
	▽ 12mm	C17	C-2	C-1	DS12
	▽ 16mm	C22	FS-3	C-12	DS16
	▽ 18mm	C6	FS-2	C-13	DS18
	▽ 24mm	C7	FS-1	C-14	DS24
	▽ 25mm	C25	-	-	-
	▽ 26mm	C26	FS	-	-
	▽ 30mm	C10	FSL	C-16	DS30
	▽ 40mm	C14	FSLX	C-17	
	▽ 50mm	C50	-	-	-
	⊽77mm	C77	LR	-	-

Circle TaperPoint

1/2 Circle Taper Point \odot

Needle	Size	SilverGlide™	Ethicon [™] (similar)	Syneture [™]	B.Braun ™ (similar)
\bigcup	⊙ 11mm	T11	-		
\bigcup	⊙ 13mm	T22	TF	CV-22	HR13
	⊙ 16mm	T16	-		
	⊙ 18mm	T28	RB-1	CV-23	HR18
	⊙ 22mm	T4	SH-1	CV-25	HR22
	⊙ 26.5mm	T5	SH /CT-1	V20/GS22	HR26
	⊙ 36mm	T8	MH /V-26	V-26	HR36
	⊙ 37mm	T37	CT-1	-	-
	⊙ 40mm	T40	-	-	-
	⊙ 65mm	Т39	TP-1	-	-

3/8Circle Taper Point \odot



Straight Cutting

Straight Cutting ∇

Needle	Size	SilverGlide™	Ethicon ™ (similar)	Syneture™	B.Braun [™] (similar)
Straight Cutting	▽ 60mm	CE60	KS	SC-2	CS60
Straight Cutting	▽ 65mm	CE65	KS	SC-2	CS60
1/2 Circle Conventional Cutting / Straight Shank (Ski) \triangle	△ 75mm	SKI75	-	-	-



Suture Material Absorption Profile

NON-ABSORBABL	E	MATERIAL DESCRIPTION	OUA	LITIES	SIMILAR TO
NYLON		Black monofilament	High strength retention		Ethilon™* Monosof™§/ Dermalon™§ Dafilon™†
FLUOROMID™		Fluorescent Yellow polyamide multifilament	Easy handling, ultra high visibili Fluorescent yellow	ty	
WHITE FLUOROMID™		White polyamide multifilament	Easy handling		
SHORTY" POLYVINYLIDENE FLUORIDEYELLOW		Bright Yellow Monofilament	High visibility Exceptional Knot Security, knot themselves, Unwetted by blood	_	
BLUE POLYPROPYLENE		Dark blue monofilament		Unwetted by blood, inert, knots flatten and lock against themselves for improved knot security	
FLUORESCENT PINK POLYPROPYLENE		Fluorescent Pink monofilament	Ultra high visibility Fluorescent pink Unwetted by blood, inert, knots flatten and lock against themselves for improved knot security		
POLYESTER		Braided polyester green silicone coated	High strength braid		Tiaron § Ethibond *
ABSORBABLE		MATERIAL DESCRIPTION	BSR PROFILE	MASS ABSORPTION	SIMILAR TO
GLYCOLIDE & ECAP		Monofilament purple	In vivo strength retention 60% at 7 days 25% at 14 days	Essentially complete between 2-4 months	Monocryl ^{™*} Caprosyn ^{™§} / Biosyn ^{™§} Monosyn ^{™†}
POLYGLYCOLICACID		Braided violet	In vivo strength retention 65% at 14 days 20% at 21 days	Essentially complete between 3-6 months	Vicryl ^{™*} Safil ^{™†} Polsorb [§]
POLYGLYCOLIC ACID FAST ABSORBING		Braided undyed	In vivo strength retention 48% at 7 days 0% at 14 days	Essentially completeat 6 weeks	Vicryl Rapide* Safil Quick†
POLYGLACTIN		Braided violet	In vivo strength retention 74% at 14 days 49% at 21 days 24% at 28 days	Essentially completeat 9 weeks	Vicryl™
POLYDIOXANONE		Monofilament purple	In vivo strength retention 80–90% at 14 days 60–80% at 28 days 50–80% at 42 days	Essentially complete between 6-9 months	PDS11™* MonoPlus™† Maxon§

^{*}Trademark of Bhicon § Trademark of Covidien † Trademark of B.Braun

Cerman Adson 1 x 2 Tissue Forcers



	and party	
Item	Product Description	UOM
54507	Cerman Adson 1 x 2 Tissue Forceps, 4.75" Length	Each



GermanAdsonBrownTissue Forceps, 4.75" USI BSI



AdsonBrown Tissue Forceps are for holding and manipulating delicate tissues. These forceps have a wide, flat thumbograsparea that is commonly serrated. The jaws are short and the tips are narrowwith 7x7 teeth.

Item	Product Description	UOM
75364G	Cerman Adson Brown ,4.75" Length	Each



GermanAllis Tissue Forceps, 6"



Cerman Allis Tissue Forceps are surgical instrument with sharp teeth, used to hold or grasp heavy tissue. Cerman Allis Tissue Forceps can cause damage, so they are often used on tissue about to be removed.

Item	Product Description	UOM
75710G	German Allis Tissue Forceps, 6"Length	Each



German Mosquito Hemostats Straight, 5"



Cerman Mosquito Hemostats can be used in various veterinary surgical and ophthalmologic procedures. They are best knownfortheir fine tips, short and fully serrated jaws. They are available with straightforcurved tips.

	- · · · · · · · · · · · · · · · · · · ·	
ltem	Product Description	UOM
75302SG	Cerman Mosquito Hernostats Straight, 5"Length	Each



German Mosquito Hemostats Curved, 5"



Cerman Mosquito Hernostats can be used in various veterinary surgical and ophthalmologic procedures. They are best knownfortheir fine tips, short and fully serrated jaws. They are available with straight or curved tips.

Ite	m	Product Description	UOM		
75	302CG	Cerman Mosquito Hemostats Curved, 5"Length	Each		



GermanKelly HemostatsStraight, 55"



Cerman Kelly Hernostats are primarily used for damping large blood vessels or manipulating heavy tissue by the veterinarian. They may also be used for soft tissue dissection. They are available curved or straight with a ratcheted finger ring handle. The tips of German Kelly Hemostats are serrated halfway up the jaws.

Item	Product Description	LOM
75306SG	GermanKelly Hemostats Straight, 5.5"Length	Each











GermanKelly Hemostats Curved 55"

USI

German Kelly Hémostats are primarily used for damping large blood vessels or manipulating heavy tissue by the veterinarian. They may also be used for soft tissue dissection. They are available curved or straight with a ratcheted finger ring handle. The tips of German Kelly Hemostats are serated halfway up the jaws.

Item	Product Description	UOM
75306CG	Cerman Kelly Hemostatis Curved, 5.5" Length	Each



Cerman Carmalt Hernostats Curved 625"



Nicknamed the "stars and stripes hemostat," German Carmalt Hemostats are characterized by the longitudinal serrations that run the length of the blade with cross-hatching at the tip. These large, crushing hemostatic forceps are a choice instrument for damping blood vessels and large tissues or ligating pedides. The textured blades, along with the scissor like ratchet handles, ensure a strong grip of the damped tissue with minimal chance of cutting the tissue cluring a procedure.

ltem	Product Description	UOM
75316CC	CermanCarmalt Hemostats Curved, 625"Length	Each



Cerman Badkhaus Towel Clamps 35"



Cerman Badkhaus Towel Clampis a perforating dampused for securing a towel or drape to the patient's skin by means of a small puncture. This product is curved with 1×1 perforating prongs.

Item	Product Description	UOM
75704G	Cerman Badkhaus Towel Clamps, 3.5" Length	Each



GermanScalpel Handle#3



GermanScalpel Handle#3 is compatible with size (10, 11, 12, 15) series disposable blades and used for sharp dissection. This handle has horizontal serrations for thumbandfore finger gripping as well as a ruler on the back side.

ltem	Product Description	UOM
75728G	CermanScalpel Handle#3	Each




Equine Eye Pack



- Plastic Sterilization Tray
- Castroviejo Eye Speculum
- Bishop Harmon1x2
- Chalazion Forceps 23 mm
- German Derf Needle Holder Tungsten Carbide
- Stevens TenotomySassors Curved

Item	Product Description	UOM
75220	USI Equine Eye Pack	Each



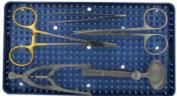

Eye Pack




- Plastic Sterilization Tray
- Castroviejo Needle Holders
- Bishop Harmon1x2
- Chalazion Forceps 20mm
- Stevens TenotomySassors Curved
- Wire Eyelid Speculum

Item	Product Description	UOM
75215	Eye Pack	Each







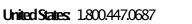



Battery-Powered Cautery Kit MMDS-VET

The Battery-Powered Cautery Kit is convenient and portable for cauterizing bleeders, removing skin tags or accomplishing many other tasks normally requiring an electro-cautery setup. At 2,200 degrees Fahrenheit, the kit is also useful in trimming beaks and nails on smaller patients and exotic animals.

Kit indudes:

Item

48426

- Fine Tip
- Fine Tip, 4"
- Fine Tip (Loop)
- 2"AA" batteries
- Carrying /storage case

Product Description

	48420	Battery-Powered Cautery Kit	Each
Replacement Tips			
	48421	Fine Tip	Each
		Fine Tip (Loop)	Each
	48423	HDTp, 018"	Each
	48424	Fine Tip, 4"	Each



WOM

Each



Fine Tip, 6"



High density thermoformed black plastic construction.

• Dimensions: 18" long x 10" wide x 12" high

- Two side compartments
- Similar to the Pandora Vet Case

Item	Product Description	UOM
01120	Vet Case, DMMMedical	Each





- Unique and specially designed outting instrument used to out and remove footwaps, ID bands and ear tags
- Designed to use on most standing cows, eliminating the need for bandage removal by restraint with a hooftrimmingtable or by lifting the bandaged leg with a rope
- Saves time and reduces the risk of injury to the bandaged patient
- Recessed cutting blade, enclosed within a guarded frame, allows user to slide the instrument under a bandaged hoof and safely out the bandage
- Durable handle grip contoured for comfortandesse of use
- Handle features a guard extending the length of the instrument providing separation of the user from the autting edge of the instrument
- Manufactured of durable plastic with the blade securely insert-molded into the handle
- Manyadditional uses including the removal of matted furon dogs and cats, efficient use in large animal surgeries or any other applications requiring safe and effective cutting
- Patent pendingglobally
- **Autodavable**











 -



Thanks for visiting Equi-Still Equine Dental Instruments!



179 Franklin Road, RD 2 Waihi 3682 NEW ZEALAND T +64 (0) 7 863 9312 www.EquiStillStocks.com