


EQUI-STILL
EQUINE STOCKS

Portable and stationary stocks...
for everyone's safety

2023 CATALOGUE



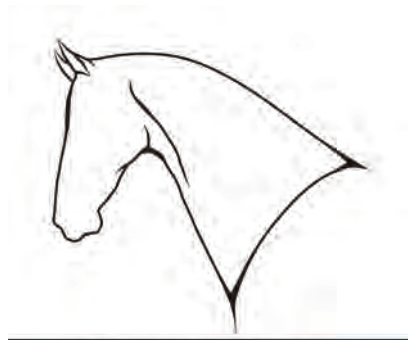

**EQUI-STILL
EQUINE STOCKS**

**PORTABLE STOCKS
TUI**

**STATIONARY STOCKS
TUI
KEA**

**OPTIONAL EXTRAS
MINI-BAR
EQUI-TABLE
EQUI-TIE
EQUI-COVER**

COMMERCIAL GUARANTEE



#1 New Zealand manufacturer of equine portable and stationary stocks

PORTABLE STOCKS

a lasting investment in your practice and your safety

- designed for safety, robustness and ease of use by a U.C. Davis-trained equine veterinarian and ICU Technician
- heavy, hot-dipped galvanized steel
- best quality closed cell foam padding with heavy duty woven polyester coated on both sides with flexible PVC
- absence of a fourth post
- quick-release rear latch
- winch, water tank, and pump. Options for power via solar charger and deep cycle battery vs. Jump-start battery, allowing you to use stocks with multiple practice vehicles without specialty wiring
- one **Equi-Table** is included, which may be mounted on one of three corners
- stocks tip up and down via winch, and may be raised and lowered even when not attached to vehicle
- patients may be restrained in forward or rearward configuration to assist in dental, reproductive, radiological, or surgical procedures
- best quality engineering: which means our stocks back and tow beautifully, barely tracking around bends ➤ our make and models are registered with the USA NHTSA (National Highway Traffic Safety Authority) ➤ stocks met voluntary NATM Compliance Certification in 2018, meaning they met or exceeded most USA states' registration requirements

MODEL TUI

- Features hydraulic brakes and a single axle
- one side of stocks is formed from square steel bars, allowing better access to patient
- opposite side is stained plywood, which may be drilled to allow rope placement to fit your smaller patients
- a rope, or Optional Extra **Equi-Tie**, holds patient in place, whether horse is placed head-first or rear-first for reproduction work



STATIONARY STOCKS

- designed for safety, robustness and ease of use by a U.C. Davis-trained equine veterinarian and ICU Technician
- heavy, hot dipped galvanized steel
- best quality closed cell foam padding with heavy duty woven polyester coated on both sides with flexible PVC
- absence of a fourth post
- quick-release rear latch, which combined with only three posts, means stocks open up fully in emergency
- stationary stocks come with separate mounting pipes to inset into floor concrete where you wish the stocks to be placed, so the stocks may later be moved
- stocks include one **Equi-Table**
- stationary stocks have no floor, winch, pump, water tank, etc.
- a lasting investment in your practice and your safety

TUI STATIONARY

- as with our portable stocks, one side of stocks is formed from square steel bars, allowing better access to patient
- a rope, or Optional Extra **Equi-Tie** holds the patient in place, whether horse is placed head-first or rear-first for reproduction work
- best quality closed cell foam padding with heavy duty woven polyester coated with flexible PVC on both sides
- absence of a fourth post • quick-release rear latch
- a lasting investment in your practice and your safety



KEA STATIONARY

- two sides of stocks are solid steel with holes for bars
- a rope, or Optional Extra **Equi-Tie**, may be used instead to hold patient in place, whether horse is placed into stocks head-first, or rear-first



OPTIONAL EXTRAS

THE MINI BAR

- mounted on the lowest side rail of the stocks to help contain your smaller patients, be they equine, camelid, bovine or anything else
- easy on, easy off




THE EQUI-TIE

- heavy, thick, soft nylon rope with brass snap
- keeps patients in place during procedures
- plywood backing of *Tui* model may be drilled at heights appropriate to your patients' sizes to run the soft rope through
- Available in 2.5 m and 5 m lengths

THE EQUI-TABLE

- removable table, one provided with stocks
- additional tables may be purchased, to your specs
- mounts on any of three different corners of stocks
- designed to tip to read ultrasound or other screens

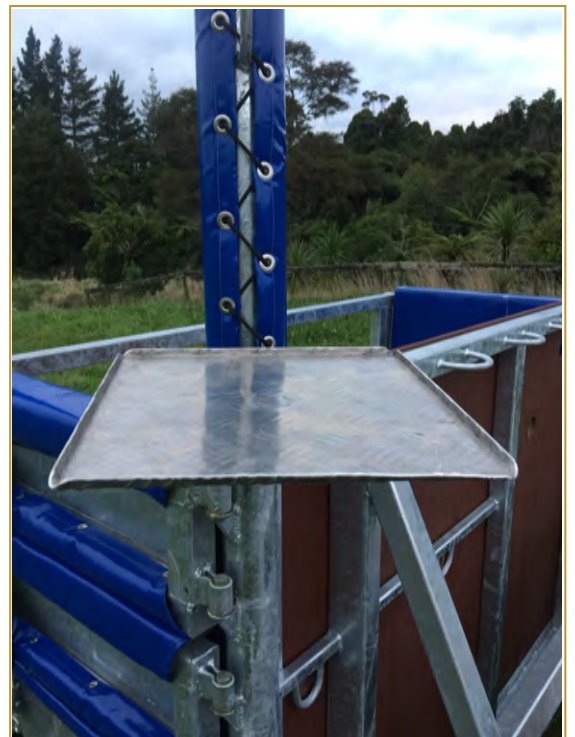

THE EQUI-COVER

- stocks cover for use while on the road or for storage
- same poly/PVC as stocks upholstery

OUR COMMERCIAL GUARANTEE

1. SCOPE

In addition to the Legal Guarantee, EQUI-STILL PORTABLE EQUINE STOCKS (EQUI-STILL) offers practitioners an EQUI-STILL Commercial Guarantee which applies to all new EQUI-STILL-branded portable and stationary stocks purchased by the practitioner from EQUI-STILL, from Aug 20, 2020.

2. DEFINITION

a) Product Commercial Guarantee

EQUI-STILL guarantees to the practitioner **against** breakage his/her EQUI-STILL-branded new portable and stationary stocks for a 12-month period following its date of purchase from EQUI-STILL. This applies to:

-Broken parts (except “wear and tear” parts - paragraph 3.

b) -Labour costs only.

b) After-Sales Service Commercial Guarantee

All repairs made by EQUI-STILL on EQUI-STILL portable and stationary stocks are guaranteed for 6 months. This guarantee covers:

-Broken parts replaced in the last 6 months (except “wear and tear” parts - paragraph 3.b)

-Labour costs for the replacement of these parts

If any repairs to EQUI-STILL Stocks have been performed without written authorisation by EQUI-STILL, this commercial guarantee is null and void.

3. IMPLEMENTATION

a) Product registration

The product commercial guarantee, the after-sales service commercial guarantee and the annual check-up will only be implemented if EQUI-STILL products are registered with EQUI-STILL within 15 days following the delivery to the practitioner from EQUI-STILL or one of its official distributors.

b) Disclaimer

The commercial guarantee applies to all EQUI-STILL Portable Equine Stocks except:

A – The following parts:

-tyres, (other than manufacturer's warranty on defects).

-jump-start battery (other than manufacturer's warranty) .

-electrical system or electrical components, if any alteration to the electrical power system or power items have been performed other than by an Official EQUI-STILL Service Center. This includes, but is not limited to, power input source, water pump, winch, lights, or any wiring.

B – If the problem is caused by an accident

C – If the damage is caused by the user, including driving at speed over irregular roads or speed bumps.

D – If the user fails to comply with the instruction manual and/or EQUI-STILL cleaning and service recommendations


E – If the VIN number of the product is removed or modified

F – If the product is modified by anyone other than an Official EQUI-STILL Service Center or if it is modified with non-EQUI-STILL parts.

G – If the product was not registered in writing by the practitioner with EQUI-STILL within 15 days after its delivery to the purchaser



Thanks for visiting
Equi-Still Stocks!

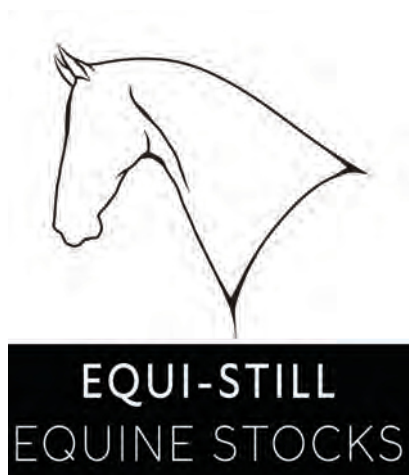

179 Franklin Road, RD 2
Waihi 3682
NEW ZEALAND
T +64 (0) 7 863 9312
www.EquiStillStocks.com